

CERTIFICATION REPORT

The certification of the catalytic activity concentration of creatine kinase in ERM®-AD455k/IFCC

European Commission
Joint Research Centre
Directorate F – Health, Consumers and Reference Materials

Contact information
Reference materials sales
Address: Retieseweg 111, 2440 Geel, Belgium
E-mail: jrc-rm-distribution@ec.europa.eu
Tel.: +32 (0)14 571 705

JRC Science Hub
<https://ec.europa.eu/jrc>

Legal Notice

This publication is a Reference Materials Report by the Joint Research Centre, the European Commission's in-house science service. It aims to provide evidence-based scientific support to the European policy-making process. The scientific output expressed does not imply a policy position of the European Commission. Neither the European Commission nor any person acting on behalf of the Commission is responsible for the use which might be made of this publication.

All images © European Union 2016

JRC 102949

EUR 28101 EN

ISBN 978-92-79-61786-7 (PDF)
ISSN 1831-9424 (online)
doi: 10.2787/310917

Luxembourg: Publications Office of the European Union, 2016
© European Union, 2016

Reproduction is authorised provided the source is acknowledged.
Printed in Belgium

Abstract

This report describes the production of ERM[®]-AD455k/IFCC, which is a material certified for the catalytic activity concentration of creatine kinase (CK). This material was produced in collaboration with the International Federation for Clinical Chemistry and Laboratory Medicine (IFCC) following ISO Guide 34:2009 and it is certified in accordance with ISO Guide 35:2006.

The starting material was a recombinant form of human CK isoenzyme MM (CK-MM) expressed in *E. coli*. It was produced, purified, filled and lyophilised by Asahi Kasei Pharma Corporation (Tokyo, Japan).

Between unit-homogeneity was quantified and stability during dispatch and storage were assessed in accordance with ISO Guide 35:2006.

The material was characterised by an interlaboratory comparison of laboratories of demonstrated competence and adhering to ISO/IEC 17025:2005. Uncertainties of the certified values were calculated in accordance with the Guide to the Expression of Uncertainty in Measurement (GUM) and include uncertainties related to possible inhomogeneity, instability and characterisation.

The material is intended for the assessment of method performance of the IFCC primary reference measurement procedure for the measurement of the catalytic activity concentration of CK at 37 °C. As with any reference material, it can be used for establishing control charts or validation studies. The CRM is available in glass vials containing lyophilised material from 1 mL of CK solution. The minimum amount of sample to be used is 12 µL after reconstitution of the whole content in one vial.

The CRM was accepted as European Reference Material (ERM[®]) after peer evaluation by the partners of the European Reference Materials consortium.

CERTIFICATION REPORT

The certification of the catalytic activity concentration of creatine kinase in ERM®-AD455k/IFCC

**Julia Kuhlmann, Brigitte Toussaint, Heinz Schimmel,
Hendrik Emons**

European Commission, Joint Research Centre
Directorate F – Health, Consumers and Reference Materials
Geel, Belgium

Disclaimer

Certain commercial equipment, instruments, and materials are identified in this paper to specify adequately the experimental procedure. In no case does such identification imply recommendation or endorsement by the European Commission, nor does it imply that the material or equipment is necessarily the best available for the purpose.

Summary

This report describes the production of ERM[®]-AD455k/IFCC, which is a material certified for the catalytic activity concentration of creatine kinase (CK). This material was produced in collaboration with the International Federation for Clinical Chemistry and Laboratory Medicine (IFCC) following ISO Guide 34:2009 [1] and it is certified in accordance with ISO Guide 35:2006 [2].

The starting material was a recombinant form of human CK isoenzyme MM (CK-MM) expressed in *E. coli*. It was produced, purified, filled and lyophilised by Asahi Kasei Pharma Corporation (Tokyo, Japan).

Between unit-homogeneity was quantified and stability during dispatch and storage were assessed in accordance with ISO Guide 35:2006 [2].

The material was characterised by an interlaboratory comparison of laboratories of demonstrated competence and adhering to ISO/IEC 17025:2005 [3].

Uncertainties of the certified values were calculated in accordance with the Guide to the Expression of Uncertainty in Measurement (GUM) [4] and include uncertainties related to possible inhomogeneity, instability and characterisation.

The material is intended for the assessment of method performance of the IFCC primary reference measurement procedure for the measurement of the catalytic activity concentration of CK at 37 °C. As with any reference material, it can be used for establishing control charts or validation studies. The CRM is available in glass vials containing lyophilised material from 1 mL of CK solution. The minimum amount of sample to be used is 12 µL after reconstitution of the whole content in one vial.

The CRM was accepted as European Reference Material (ERM[®]) after peer evaluation by the partners of the European Reference Materials consortium.

The following value was assigned:

	Certified value ²⁾	Uncertainty ³⁾
Catalytic activity concentration ¹⁾	314 U/L 5.2 µkat/L	6 U/L 0.1 µkat/L
¹⁾ Catalytic activity concentration of creatine kinase isoenzyme MM (CK-MM) in the reconstituted material, as obtained by the IFCC primary reference measurement procedure for the measurement of catalytic activity concentration of creatine kinase at 37 °C.		
²⁾ Certified values are values that fulfil the highest standards of accuracy and represent the unweighted mean value of the means of accepted sets of data, each set being obtained in a different laboratory. The certified value and its uncertainty are traceable to the International System of Units (SI). Values were converted from U/L into µkat/L by multiplication with the factor $f = 0.01667$.		
³⁾ The uncertainty is the expanded uncertainty of the certified value with a coverage factor $k = 2$ corresponding to a level of confidence of about 95 % estimated in accordance with ISO/IEC Guide 98-3, Guide to the Expression of Uncertainty in Measurement (GUM:1995), ISO, 2008.		

Table of contents

Summary.....	1
Table of contents.....	2
Glossary.....	4
1 Introduction.....	7
1.1 Background	7
1.2 Choice of the material	7
1.3 Design of the CRM project	8
2 Participants	8
2.1 Project management and evaluation	8
2.2 Processing.....	8
2.3 Homogeneity and stability study	8
2.4 Characterisation.....	8
3 Material processing and process control	9
3.1 Origin of the starting material and processing	9
3.2 Process control	9
4 Homogeneity	9
4.1 Between-unit homogeneity.....	10
4.2 Within-unit homogeneity and minimum sample intake.....	11
5 Stability.....	11
5.1 Short-term stability study	12
5.2 Long-term stability study	12
5.3 Estimation of uncertainties	13
6 Characterisation	14
6.1 Selection of participants.....	14
6.2 Study setup.....	14
6.3 Methods used	14
6.4 Evaluation of results	14
6.4.1 Technical evaluation	14
6.4.2 Statistical evaluation	15
7 Value Assignment.....	15
7.1 Certified values and their uncertainties	16
8 Metrological traceability and commutability.....	16
8.1 Metrological traceability	16
8.2 Commutability	17
9 Instructions for use	17
9.1 Safety information.....	17
9.2 Storage conditions	17
9.3 Reconstitution	18
9.4 Minimum sample intake	18
9.5 Use of the certified value	18
10 Acknowledgments	20
11 References	21
Annexes	23

Glossary

ANOVA	Analysis of variance
CK-MM	Creatine kinase MM isoenzyme
CRM	Certified reference material
C-RSE	IFCC Committee for Reference Systems of Enzymes
EC	European Commission
<i>E. coli</i>	<i>Escherichia coli</i>
ERM [®]	Trademark of European Reference Materials
EU	European Union
GUM	Guide to the Expression of Uncertainty in Measurements [ISO/IEC Guide 98-3:2008]
IFCC	International Federation of Clinical Chemistry and Laboratory Medicine
IRMM	Institute for Reference Materials and Measurements of the JRC
ISO	International Organization for Standardization
k	Coverage factor
kat/L	Katal per litre
m	Mass
MS_{between}	Mean of squares between-unit from an ANOVA
MS_{within}	Mean of squares within-unit from an ANOVA
n	Number of replicates per unit
p	Number of datasets
rel	Index denoting relative figures (uncertainties etc.)
s	Standard deviation
s_{bb}	Between-unit standard deviation; an additional index "rel" is added when appropriate
s_{between}	Standard deviation between groups as obtained from ANOVA; an additional index "rel" is added as appropriate
SI	International System of Units
s_{rel}	Relative standard deviation of all results of the stability study
s_{within}	Standard deviation within groups as obtained from ANOVA; an additional index "rel" is added as appropriate
s_{wb}	Within-unit standard deviation
\bar{t}	Mean of all t_i
t_i	Time point for each replicate
t_{sl}	Proposed shelf life
t_{tt}	Proposed transport time
u	Standard uncertainty

U	Expanded uncertainty
U/L	Units per litre
u_{bb}^*	Standard uncertainty related to a maximum between-unit inhomogeneity that could be hidden by method repeatability/intermediate precision select as appropriate; an additional index "rel" is added as appropriate
u_{bb}	Standard uncertainty related to a possible between-unit inhomogeneity; an additional index "rel" is added as appropriate
u_{char}	Standard uncertainty of the material characterisation; an additional index "rel" is added as appropriate
u_{CRM}	Combined standard uncertainty of the certified value; an additional index "rel" is added as appropriate
U_{CRM}	Expanded uncertainty of the certified value; an additional index "rel" is added as appropriate
u_{Δ}	Combined standard uncertainty of measurement result and certified value
U_{Δ}	Expanded standard uncertainty of measurement result and certified value
u_{lts}	Standard uncertainty of the long-term stability; an additional index "rel" is added as appropriate
u_{meas}	Standard measurement uncertainty
u_{sts}	Standard uncertainty of the short-term stability; an additional index "rel" is added as appropriate
v	Volume
\bar{x}	Arithmetic mean
Δ_{meas}	Absolute difference between mean measured value and the certified value
$\nu_{MS_{within}}$	Degrees of freedom of MS_{within}

1 Introduction

1.1 Background

Creatine kinase (CK) is an enzyme that is generally associated with adenosine triphosphate regeneration in contractile or transport systems. Its predominant physiologic function occurs in muscle cells, where it is involved in the storage and utilisation of high-energy creatine phosphate. It is widely distributed in tissues, with highest activities found in skeletal muscle, heart muscle, and brain tissue. Creatine kinase is a dimeric molecule consisting of two elongated polypeptide subunits termed B (brain type) and M (muscle type) with three known cytosolic isoenzymes (MM, MB and BB). The cytosolic isoenzymes can also be named on the basis of their electrophoretic mobility towards the anode [5] and are thus called: creatine kinase-1 (CK-BB), creatine kinase-2 (CK-MB) and creatine kinase-3 (CK-MM). A fourth isoenzyme occurs in mitochondria of several tissues and is composed of a third type of subunit [6]. Creatine kinase gene expression is tissue specific as well as developmentally and hormonally regulated [7]. Measurement of this enzyme in serum is requested to the clinical laboratory most frequently for the diagnosis of acute myocardial infarction.

The catalytic activity of an enzyme is a property that is measured by the rate of a specified chemical reaction under certain experimental conditions. The measurement of this property is very important in clinical chemistry, but the standardisation of catalytic activity measurements is challenging as a number of parameters influence the enzyme activity (e.g. temperature, pH, substrate nature and concentration, activators, inhibitors). Therefore, the measurement results are heavily dependent on the measurement procedure used to attain them. This led to the development of universally recognised measurement procedures for enzymes commonly measured in clinical chemistry, such as the IFCC primary reference measurement procedure for the measurement of catalytic activity concentrations of enzymes at 37 °C [8].

The EU directive on *in vitro* diagnostic medical devices (Directive 98/79/EC) requires traceability of the assigned values of calibrants and control materials to reference measurement procedures and/or reference materials of higher order.

In collaboration with the IFCC Committee for Reference Systems of Enzymes (C-RSE), the Institute for Reference Materials and Measurements (IRMM) developed a CRM certified for the catalytic activity concentration of CK. This material, ERM-AD455k/IFCC, is intended to be used as a quality control material for the IFCC primary reference measurement procedure for CK at 37 °C [9]. The homogeneity and the stability of ERM-AD455k/IFCC were demonstrated and the certified value was assigned using the IFCC reference measurement procedure at 37 °C in an interlaboratory comparison of expert laboratories.

1.2 Choice of the material

A recombinant CK-MM material from Asahi Kasei Pharma Corporation was selected as starting material based on the outcome of a commutability study carried out by IRMM. The purity was guaranteed by the provider as containing no other enzyme with a relative catalytic activity concentration of more than 1.0 % of the total catalytic activity concentration. The material was solubilised in a buffer, frozen and lyophilised to improve long-term stability. The aim of the production process was to obtain a material that once reconstituted with 1.0 mL of distilled/deionised water would have a catalytic activity concentration of about 300 U/L. The selection of the material and the catalytic activity concentration was based on a recommendation received from the IFCC C-RSE in 2011.

1.3 Design of the CRM project

A commutability study including nine routine measurement procedures was completed to select the most appropriate starting material for the production of ERM-AD455k/IFCC.

The material was certified by interlaboratory comparison of data from expert laboratories using the IFCC primary reference measurement procedure for the measurement of the catalytic activity concentration of CK at 37 °C [9]. The performances of the laboratories were assessed using two serum based control materials.

The homogeneity and stability of the material were assessed using a UniCel® DxC 800 Synchron Clinical System with CK reagent cartridges (Beckman Coulter, Inc., Clare, IE). This test kit showed a low average relative standard deviation during the commutability study. Statistical analysis of the data was done using the software SoftCRM (version 2.0.21) [10].

2 Participants

2.1 Project management and evaluation

European Commission, Joint Research Centre, Institute for Reference Materials and Measurements (IRMM), Geel, BE
(accredited to ISO Guide 34 for production of certified reference materials, BELAC No. 268-RM)

2.2 Processing

Asahi Kasei Pharma Corporation, Tokyo, JP

2.3 Homogeneity and stability study

European Commission, Joint Research Centre, Institute for Reference Materials and Measurements (IRMM), Geel, BE
(accredited to ISO Guide 34 for production of certified reference materials, BELAC No. 268-RM)

Beckman Coulter, Inc., Clare, IE

2.4 Characterisation

The laboratories are listed below in alphabetic order. This order does not necessarily correspond to the ranking of the laboratories L01 to L10 described in the Tables of this report.

*Affiliated Hospital of Nantong University, Reference Laboratory, Nantong, CN
(measurements under the scope of ISO/IEC 17025:2005 and ISO 15195:2003 accreditation, CNAS No. L6260)

*Beijing Aerospace General Hospital, Reference Laboratory, Beijing, CN
(measurements under the scope of ISO/IEC 17025:2005 and ISO 15195 accreditation, CNAS No. L5536)

Biosystems, S.A., Barcelona, ES

*INSTAND e.V., Düsseldorf, DE
(measurements under the scope of ISO/IEC 17025:2005 accreditation, DAkkS No. D-K-15027-01-00)

LabWest HagaZiekenhuis, Klinisch Chemisch en Hematologisch Laboratorium, Den Haag, NL

Servizio Di Medicina Di Laboratorio Diagnostica e Ricerca, Ospedale San Raffaele, Milano, IT

*Sichuan Maker Biotechnology Co., Ltd., Reference System Department, Chengdu, CN
(measurements under the scope of ISO/IEC 17025:2005 and ISO 15195:2003 accreditation, CNAS No. L6172)

*Stiftung für Pathobiochemie und Molekulare Diagnostik, Referenzinstitut für Bioanalytik, Kalibrierlaboratorium II, Medizinische Hochschule Hannover, Institut für Klinische Chemie, Hannover, DE
(measurements under the scope of ISO/IEC 17025:2005 and ISO 15195:2004 accreditation, DAkkS No. D-K-15117-02-00)

*Universitat Autònoma de Barcelona, Departament de Bioquímica i Biologia Molecular, Unitat de Bioquímica de Medicina, Barcelona, ES
(measurements under the scope of ISO/IEC 17025:2005 and ISO 15195:2004 accreditation, ENAC No. 195/LC10.141)

*Università degli Studi di Milano, Laboratorio Analisi Chimico-Cliniche, Centro Interdipartimentale di Ricerca sulla Riferibilità Metrologica in Medicina di Laboratorio, Milano, IT
(measurements under the scope of ISO/IEC 17025 and ISO 15195:2004 accreditation, ACCREDIA No. 217/01)

*Listed in the JCTLM database for reference measurement services

3 Material processing and process control

3.1 Origin of the starting material and processing

The starting material was produced and processed by Asahi Kasei Pharma Corporation. It was a recombinant form of human CK-MM expressed in *E. coli*. The enzyme was purified by hydrophobic and anion-exchange chromatography and then dissolved in a buffer (pH 7.5) containing among others bovine serum albumin and polysaccharides. The liquid preparation was filled in glass vials, frozen and lyophilised. The vials were labelled in the same order as the filling sequence using the labels provided by IRMM.

3.2 Process control

The final starting material was prepared and tested for activity, stability and homogeneity by Asahi Kasei Pharma Corporation. This batch was shown to be homogenous and stable for at least one year when stored at -20 °C.

4 Homogeneity

A key requirement for any reference material aliquotted into units is equivalence between those units. In this respect, it is relevant whether the variation between units is significant compared to the uncertainty of the certified value, but it is not relevant if this variation between units is significant compared to the analytical variation. Consequently, ISO Guide 34 [1] requires reference material producers to quantify the between-unit variation. This aspect is covered in between-unit homogeneity studies.

The within-unit inhomogeneity does not influence the uncertainty of the certified value when the minimum sample intake is respected, but determines the minimum size of an aliquot that is representative for the whole unit.

4.1 Between-unit homogeneity

The between-unit homogeneity was evaluated to ensure that the certified value of the CRM is valid for all vials of the material, within the stated uncertainty.

The number of vials selected corresponds to approximately the cube root of the total number of vials produced. The 20 vials were selected using a random stratified sampling scheme covering the whole batch for the between-unit homogeneity test. For this, the batch was divided into 20 groups (with a similar number of vials) and one vial was selected randomly from each group. After reconstitution, three samples were taken from each selected vial, and analysed with a UniCel® DxC 800 Synchron Clinical System using CK reagent cartridges. The measurements were performed under repeatability conditions, and in a randomised manner to be able to separate a potential trend in the analytical sequence from a trend in the filling sequence.

Regression analyses were performed to evaluate potential trends in the analytical sequence as well as trends in the filling sequence. No trends in the filling sequence or the analytical sequence were observed at a 95 % confidence level.

The dataset was assessed for consistency using Grubbs outlier tests at a confidence level of 99 % on the individual results and on the unit means. No outlying individual results and outlying unit means were detected. All the data were retained for statistical analysis.

Quantification of between-unit inhomogeneity was undertaken by analysis of variance (ANOVA), which separates the between-unit variation (s_{bb}) from the within-unit variation (s_{wb}). The latter is equivalent to the method repeatability if the individual samples were representative for the whole vial.

Evaluation by ANOVA requires mean values per vial, which follow at least a unimodal distribution and results for each vial that follow unimodal distributions with approximately the same standard deviations. The distribution of the mean values per vial was visually tested using histograms and normal probability plots.

It should be noted that $s_{bb,rel}$ and $s_{wb,rel}$ are estimates of the standard deviations and are therefore subject to random fluctuations. Therefore, the mean square between groups ($MS_{between}$) can be smaller than the mean squares within groups (MS_{within}), resulting in negative arguments under the square root used for the estimation of the between-unit variation, whereas the true variation cannot be lower than zero. In this case, u_{bb}^* , the maximum inhomogeneity that could be hidden by method repeatability, was calculated as described by Linsinger *et al.* [11]. u_{bb}^* is comparable to the limit of detection of an analytical method. It describes the maximum inhomogeneity that might be hidden in the frame of the given study setup.

Method repeatability ($s_{wb,rel}$), between-unit standard deviation ($s_{bb,rel}$) and $u_{bb,rel}^*$ were calculated as:

$$s_{wb,rel} = \frac{\sqrt{MS_{within}}}{\bar{x}} \quad \text{Equation 1}$$

$$s_{bb,rel} = \frac{\sqrt{\frac{MS_{between} - MS_{within}}{n}}}{\bar{x}} \quad \text{Equation 2}$$

$$u_{bb,rel}^* = \frac{\sqrt{\frac{MS_{within}}{n}} \sqrt[4]{\frac{2}{v_{MS_{within}}}}}{\bar{x}} \quad \text{Equation 3}$$

MS_{within} mean of squares within-unit from an ANOVA

$MS_{between}$ mean of squares between-unit from an ANOVA

\bar{x}	mean of all results of the homogeneity study
n	mean number of replicates per unit
$\nu_{MS_{within}}$	degrees of freedom of MS_{within}

The results of the evaluation of the between-unit variation are summarised in Table 1. The resulting values from the above equations were converted into relative uncertainties.

Table 1: Results of the homogeneity study

ERM-AD455k/IFCC	$s_{wb,rel}$ [%]	$s_{bb,rel}$ [%]	$u_{bb,rel}^*$ [%]
Catalytic activity concentration of CK	0.5	0.5	0.1

The homogeneity study showed no outlying unit means or trends in the filling sequence. Therefore the between-unit standard deviation can be used as estimate of u_{bb} . As u_{bb}^* sets the limits of the study to detect inhomogeneity, the larger value of s_{bb} and u_{bb}^* is adopted as uncertainty contribution to account for potential inhomogeneity.

4.2 Within-unit homogeneity and minimum sample intake

The within-unit homogeneity is closely correlated to the minimum sample intake. The minimum sample intake is the minimum amount of sample that is representative for the whole unit and thus should be used in an analysis. Using sample sizes equal or above the minimum sample intake guarantees the certified value within its stated uncertainty.

Homogeneity and stability experiments were performed using a 12 μ L sample intake. This sample intake gives acceptable repeatability, demonstrating that the within-unit inhomogeneity no longer contributes to analytical variation at this sample intake.

5 Stability

Stability testing is necessary to establish the conditions for storage (long-term stability) as well as the conditions for dispatch of the materials to the customers (short-term stability).

Time, temperature, light (including ultraviolet radiation) and water content were regarded as the most relevant influences on the stability of the materials. The influence of ultraviolet or visible light was minimised by storing the material in the dark and dispatched in boxes, thus removing any possibility of degradation by light. The water content was reduced by freeze-drying to obtain a stable material. Therefore, only the influences of time and temperature needed to be investigated.

During transport, especially in summer time, temperatures up to 60 °C can be reached and stability under these conditions must be demonstrated, if the samples are to be transported without any additional cooling.

The stability studies were carried out using an isochronous design [12]. In this approach, samples were stored for a particular length of time at different temperature conditions. Afterwards, the samples were moved to conditions where further degradation can be assumed to be negligible (reference conditions). At the end of the isochronous storage, the samples were analysed simultaneously under repeatability conditions which greatly improves the sensitivity of the stability tests.

5.1 Short-term stability study

For the short-term stability study, samples were stored at -20 °C, 4 °C and 60 °C for 0, 1, 2 and 4 weeks (at each temperature). The reference temperature was set to -70 °C. Two vials per storage time were selected using a random stratified sampling scheme. After reconstitution, three samples were measured from each vial with a UniCel DxC 800 Synchron Clinical System using CK reagent cartridges. The measurements were performed under repeatability conditions, and a randomised sequence was used to differentiate any potential analytical drift from a trend over storage time.

The data were evaluated individually for each temperature. The results were screened for outliers using the single and double Grubbs test on a confidence level of 99 %. One outlying individual result was found (Table2). As no technical reason for the outlier could be found all data were retained for statistical analysis.

In addition, the data were evaluated against storage time, and regression lines of catalytic activity concentration versus time were calculated, to test for potential increases/decrease of the measurand due to shipping conditions. The slopes of the regression lines were tested for statistical significance.

The results of the measurements are shown in Annex B. The results of the statistical evaluation of the short-term stability are summarised in Table 2.

Table 2: Results of the short-term stability test

ERM-AD455k/IFCC	Number of individual outlying results*			Significance of the trend**		
	-20 °C	4 °C	60 °C	-20 °C	4 °C	60 °C
Catalytic activity concentration of CK	none	1-statistical reason (retained)	none	no	no	yes

* 99 % confidence level

** 95 % confidence level

Only one statistical outlier was detected on a 99 % confidence level and retained for the estimation of u_{STS} . The only trend statistically significant on a 95 % confidence level was found for storage at 60 °C.

Although the material is stable at 4 °C for 4 weeks, it shall be shipped on dry ice to avoid the introduction of a temperature shift (influence on catalytic activity concentration not evaluated) since the long-term storage temperature is -20 °C.

5.2 Long-term stability study

For the long-term stability study, samples were stored at -20 °C and 4 °C for 0, 4, 8 and 12 months (at each temperature). The reference temperature was set to -70 °C. Two samples per storage time were selected using a random stratified sampling scheme. After reconstitution, three samples were measured from each vial with a UniCel DxC 800 Synchron Clinical System using CK reagent cartridges. The measurements were performed under repeatability conditions, in a random sequence to be able to separate any potential analytical drift from a trend over storage time.

The long-term stability data were evaluated individually for each temperature. The results were screened for outliers using the single and double Grubbs test at a confidence level of 99 %. No outlying individual results were found and all data were retained for statistical analysis.

In addition, the data were plotted against storage time and linear regression lines of catalytic activity concentration versus time were calculated. The slopes of the regression lines were tested for statistical significance (loss/increase due to storage). A significant trend was detected at a 95 % confidence level for storage at 4 °C.

The results of the long-term stability measurements are shown in Annex C.

The material shall be stored at -20 °C.

5.3 Estimation of uncertainties

Since there are always variations on measurement results, no study can entirely rule out degradation of materials, even in the absence of statistically significant trends. It is therefore necessary to quantify the potential degradation. This degradation could be hidden in the frame of the study set up, i.e. to estimate the uncertainty of stability. This means that, even under ideal conditions, the outcome of a stability study can only be that there is no detectable degradation within an uncertainty to be estimated.

The uncertainties of stability during dispatch and storage were estimated, as described in [13]. In this approach, the uncertainty of the linear regression line with a slope of zero was calculated. The uncertainty contributions u_{sts} and u_{lts} were calculated as the product of the chosen transport time/shelf life and the uncertainty of the regression lines as:

$$u_{sts,rel} = \frac{s_{rel}}{\sqrt{\sum (t_i - \bar{t})^2}} \cdot t_{tt} \quad \text{Equation 4}$$

$$u_{lts,rel} = \frac{s_{rel}}{\sqrt{\sum (t_i - \bar{t})^2}} \cdot t_{sl} \quad \text{Equation 5}$$

s_{rel} relative standard deviation of all results of the stability study

t_i time elapsed at time point i

\bar{t} mean of all t_i

t_{tt} chosen transport time (1 week at -20 °C)

t_{sl} chosen shelf life (12 months at -20 °C)

The following uncertainties were estimated:

- $u_{sts,rel}$, the uncertainty of degradation during dispatch. This was estimated from the short-term stability studies at -20 °C. The uncertainty describes the possible change during a dispatch at -20 °C lasting for one week.
- $u_{lts,rel}$, the stability during storage. This uncertainty contribution was estimated from the long-term stability studies at -20 °C. The uncertainty contribution describes the possible degradation during 12 months storage at -20 °C.

The results of these evaluations are summarised in Table 3.

Table 3: Uncertainties of stability during dispatch and storage. $u_{sts,rel}$ was calculated for a temperature of -20 °C and 1 week; $u_{lts,rel}$ was calculated for a storage temperature of -20 °C and 12 months

ERM-AD455k/IFCC	$u_{sts,rel}$ [%]	$u_{lts,rel}$ [%]
Catalytic activity concentration of CK	0.1	0.5

The material showed significant degradation at 60 °C, but no significant degradation was observed for transport below 4 °C. Transport on dry ice is therefore necessary.

After the certification study, the material will be included in the IRMM's regular stability monitoring programme, to control its further stability.

6 Characterisation

The material characterisation is the process of determining the property value of a reference material.

This was based on an interlaboratory comparison of data from expert laboratories, i.e. the catalytic activity concentration of CK in the material was determined in different laboratories. All participants used the IFCC primary reference measurement procedures for the measurement of catalytic activity concentration of CK at 37 °C [9] for the measurements. Using an interlaboratory comparison aims at randomisation of laboratory bias, which reduces the combined uncertainty.

6.1 Selection of participants

Ten laboratories were selected based on criteria that comprised both technical competence and quality management aspects. Each participant was required to operate a quality system and to deliver documented evidence of its laboratory proficiency in the field of catalytic activity concentration measurements of enzymes using the IFCC primary reference measurement procedure for the measurement of catalytic activity concentrations at 37 °C. Having a formal accreditation was not mandatory, but meeting the requirements of ISO/IEC 17025 was obligatory. Where measurements are covered by the scope of accreditation, the accreditation number is stated in the list of participants (Section 2).

6.2 Study setup

Each laboratory received six vials of ERM-AD455k/IFCC and was requested to provide six independent results, one per vial. The vials for material characterisation were selected using a random stratified sampling scheme and covered the whole batch. The sample preparations and measurements had to be spread over at least two days to ensure intermediate precision conditions. Each participant was also required to analyse two blinded quality control samples.

Laboratories were also requested to give estimations of the expanded uncertainties of the mean value of the six results. All GUM approaches were regarded as equally valid procedures.

6.3 Method used

All laboratories used the IFCC primary reference measurement procedure for the measurement of the catalytic activity concentration of CK at 37 °C [9,14].

6.4 Evaluation of results

The characterisation study resulted in eight accepted datasets. All accepted individual results of the participants are displayed in tabular and graphical form in Annex D.

6.4.1 Technical evaluation

The following criteria were considered during the evaluation:

- compliance with the analysis protocol: sample preparations and measurements performed on two days
- performance in measuring the quality control samples

Based on the above criteria, the dataset from one laboratory (L05) was not included in the evaluation based on their performance in measuring the quality control samples.

6.4.2 Statistical evaluation

The datasets accepted based on technical reasons were tested for normality of dataset means using kurtosis/skewness tests and normal probability plots and were tested for outlying means using the Grubbs test and using the Cochran test for outlying standard deviations, (both at a 99 % confidence level). Standard deviations within (s_{within}) and between (s_{between}) laboratories were calculated using one-way ANOVA. The results of these evaluations are shown in Table 4.

Table 4: Statistical evaluation of the technically accepted datasets for ERM-AD455k/IFCC. p : number of technically valid datasets

ERM-AD455k/IFCC	p	Outliers		Normally distributed	Statistical parameters			
		Means	Variances		Mean [U/L]	s [U/L]	s_{between} [U/L]	s_{within} [U/L]
Catalytic activity concentration of CK	9	none	none	yes	313.48	6.32	5.87	2.48

The laboratory means follow normal distributions. None of the data contains outlying means and variances. The datasets are therefore consistent and the mean of laboratory means is a good estimate of the true value. Standard deviations between laboratories are considerably larger than the standard deviation within laboratories, showing that confidence intervals of replicate measurements are unsuitable as estimate of measurement uncertainty.

The uncertainty related to the characterisation is estimated as the standard error of the mean of laboratory means (Table 5).

Table 5: Uncertainty of characterisation for ERM-AD455k/IFCC

ERM-AD455k/IFCC	p	Mean [U/L]	s [U/L]	$u_{\text{char,rel}}$ [%]
Catalytic activity concentration of CK	9	314	6	0.6

7 Value Assignment

Based on the results of the characterisation study a certified value for the catalytic activity concentration of CK was assigned to ERM-AD455k/IFCC.

Certified values are values that fulfil the highest standards of accuracy. Procedures at IRMM require generally pooling of not less than six datasets to assign certified values. Full uncertainty budgets in accordance with the 'Guide to the Expression of Uncertainty in Measurement' [4] were established.

7.1 Certified values and their uncertainties

The unweighted mean of the means of the accepted datasets as shown in Table was assigned as certified value.

The assigned uncertainty consists of uncertainties relating to characterisation, u_{char} (Section 6), potential between-unit inhomogeneity, u_{bb} (Section 4.1), and potential degradation during transport, u_{sts} , and long-term storage, u_{its} (Section 5). These different contributions were combined to estimate the relative expanded uncertainty of the certified value ($U_{\text{CRM,rel}}$) with a coverage factor k given as:

$$U_{\text{CRM,rel}} = k \cdot \sqrt{u_{\text{bb,rel}}^2 + u_{\text{sts,rel}}^2 + u_{\text{its,rel}}^2 + u_{\text{char,rel}}^2} \quad \text{Equation 6}$$

- u_{char} was estimated as described in Section 6
- u_{bb} was estimated as described in Section 4.1
- u_{sts} and u_{its} were estimated as described in section 5.3

A coverage factor k of 2 was applied, to obtain the expanded uncertainty. The certified value and its uncertainties are summarised in Table 6.

Table 6: Certified value and its uncertainties for ERM-AD455k/IFCC

ERM-DA455k/IFCC	Certified value [U/L]	$u_{\text{char,rel}}$ [%]	$u_{\text{bb,rel}}$ [%]	$u_{\text{sts,rel}}$ [%]	$u_{\text{its,rel}}$ [%]	$U_{\text{CRM}}^{1)}$ [U/L]
Catalytic activity concentration of CK	314	0.6	0.5	0.1	0.5	6

¹⁾ Expanded ($k = 2$) and rounded uncertainty.

The International Union of Pure and Applied Chemistry and the International Union of Biochemistry recommended that enzyme concentration is expressed in terms of katal per liter (kat/L) [14]. This name and symbol were approved by the General Conference on Weights and Measures and is consistent with the International System of Units (SI) as $\text{kat} = \text{mol/s}$. However, different units have been introduced in the past. Therefore, the Commission on Enzymes of the International Union of Biochemistry proposed the term international unit (U) as the quantity of enzyme that catalyses the reaction of 1 μmol of substrate per minute. Catalytic activity concentration is then to be expressed in terms of U/L [15, 16]. Enzyme units (U) are still more commonly used than the katal, especially in biochemistry. Therefore, the certified value for ERM-AD455k/IFCC is expressed both in $\mu\text{kat/L}$ and in U/L in this certification report as well as on the certificate. The catalytic activity concentration in $\mu\text{kat/L}$ can easily be converted to U/L by multiplying with the factor $f = 60$.

1 $\mu\text{kat/L} = 60 \text{ U/L}$

1 U = $10^{-6} \text{ mol/60 s} = 16.7 \times 10^{-9} \text{ mol/s}$

8 Metrological traceability and commutability

8.1 Metrological traceability

Identity

The catalytic activity concentration of CK is a method-defined measurand and can only be obtained by following the procedure specified in the IFCC primary reference measurement procedure at 37 °C [9]. Adherence to this procedure was confirmed by agreement of the laboratories' results with the assigned value for the samples that were used as quality control samples and by comparison among laboratory results. The measurand is therefore operationally defined by *method*.

Quantity value

Traceability of the obtained results is based on the traceability of all relevant input factors. Instruments in individual laboratories were verified and calibrated with measurement standards ensuring traceability to the SI. Consistency in the interlaboratory comparison demonstrates that all relevant input factors were covered. As the assigned value is a combination of agreeing results individually traceable to the SI, the assigned quantity value itself is traceable to the SI.

8.2 Commutability

Many measurement procedures include one or more steps which select specific (or specific groups of) analytes from the sample for the subsequent whole measurement process. Often the complete identity of these 'intermediate analytes' is not fully known or taken into account. Therefore, it is difficult to mimic all analytically relevant properties of real samples within a CRM. The degree of equivalence in the analytical behaviour of real samples and a CRM with respect to various measurement procedures (methods) is summarised in a concept called 'commutability of a reference material'. There are various definitions that define this concept. For instance, the Clinical and Laboratory Standards Institute Guideline C53-A [17] recommends the use of the following definition for the term *commutability*:

"The equivalence of the mathematical relationships among the results of different measurement procedures for an RM and for representative samples of the type intended to be measured."

The commutability of a CRM defines its fitness for use and is therefore a crucial characteristic when applying different measurement methods. When the commutability of a CRM is not established, the results from routinely used measurement procedures cannot be legitimately compared with the certified value to determine whether a bias does not exist in calibration, nor can the CRM be used as a calibrant. For instance, CRMs intended to be used to establish or verify metrological traceability of routine clinical measurement procedures must be commutable for the routine clinical measurement procedures for which they are intended to be used.

A commutability study was carried out on a trial batch of the starting material for ERM-AD455k/IFCC in collaboration with the IFCC C-RSE. The results were convincing enough to process the final batch of the material and to certify its CK catalytic activity concentration. However, if ERM-AD455k/IFCC would be used to calibrate routine procedures, extended commutability studies should be performed with these procedures.

9 Instructions for use

9.1 Safety information

The usual laboratory safety measures apply.

The material is for *in vitro* use only.

9.2 Storage conditions

Unopened vials of the material should be stored at $(-20 \pm 5) ^\circ\text{C}$ in the dark. After reconstitution, the material must be kept cold ($2-8 ^\circ\text{C}$) and must be used within four hours.

Please note that the European Commission cannot be held responsible for changes that happen during storage of the material at the customer's premises, especially for opened vials.

9.3 Reconstitution

To prepare ERM-AD455k/IFCC for use, the lyophilised material shall be reconstituted according to the following procedure:

- 1) Remove vial from freezer and let equilibrate to room temperatures (20-25 °C).
- 2) Tap the vertically positioned vial gently to ensure that the lyophilised material is at the bottom of the vial.
- 3) Carefully open vial, avoiding the loss of lyophilised material.
- 4) Weigh the vial with its content to the nearest 0.1 mg.
- 5) Reconstitute with (1.00 ± 0.01) mL distilled water (20-22 °C) slowly added to the sides of the vial.
- 6) Weigh the vial after adding the water and record the weight.
- 7) Carefully close the vial.
- 8) Allow to stand at room temperature for ten minutes.
- 9) Slowly stir up the vial to dissolve the lyophilised material completely.
- 10) Calculate the volume of water at 20 °C from the mass of water added taking into account the temperature dependent density.
- 11) Keep the vial cold (2-8 °C) until use.
- 12) The catalytic activity concentration of CK must be measured within four hours following the reconstitution*.

*The activity is not guaranteed after four hours after reconstitution.

9.4 Minimum sample intake

The minimum sample intake representative for the catalytic activity concentration of CK in ERM-AD455k/IFCC is 12 µL as this was the sample intake for the homogeneity study.

9.5 Use of the certified value

The main purpose of the material is to control the performance of the IFCC primary reference procedure for the measurement of catalytic activity concentration of CK at 37 °C [9]. As any reference material, it can be used for establishing control charts or validation studies.

Use as a calibrant

It is not recommended to use this material as calibrant. If used nevertheless, the uncertainty of the certified value shall be taken into account in the estimation of the measurement uncertainty. When the material is used as a calibrant in a routine measurement procedure the commutability should be verified for the assay concerned.

Comparing an analytical result with the certified value

A result is unbiased if the combined standard uncertainty of measurement and certified value covers the difference between the certified value and the measurement result (see also ERM Application Note 1, www.erm-crm.org [18]).

When assessing the method performance, the measured values of the CRMs are compared with the certified values. The procedure is summarised here:

- Calculate the absolute difference between mean measured value and the certified value (Δ_{meas}).
- Combine the measurement uncertainty (u_{meas}) with the uncertainty of the certified value (u_{CRM}): $u_{\Delta} = \sqrt{u_{\text{meas}}^2 + u_{\text{CRM}}^2}$
- Calculate the expanded uncertainty (U_{Δ}) from the combined uncertainty (u_{Δ}) using an appropriate coverage factor, corresponding to a level of confidence of approximately 95 %.

- If $\Delta_{\text{meas}} \leq U_{\Delta}$ then no significant difference exists between the measurement result and the certified value, at a confidence level of approximately 95 %.

Use in quality control charts

The material can be used for quality control charts. Using CRMs for quality control charts has the added value that a trueness assessment is built into the chart.

10 Acknowledgments

The authors would like to acknowledge the support received from Hanne Leys from IRMM relating to the processing of this CRM and from Maria Concepcion Contreras concerning the set-up of the required isochronous studies.

Furthermore, the authors would like to thank Virginie Tregoeat (IRMM) for reviewing the certification report, as well as the experts of the Certification Advisory Panel "Biological Macromolecules and Biological/Biochemical Parameters", Hez Hird (Food and Environmental Research Agency, UK), Martin Wagner (University for Veterinary Medicine Vienna, AT), Lothar Siekmann (University Bonn, DE) for their constructive comments. The authors are also grateful for the active involvement of the participating laboratories in the characterisation study and have appreciated the thoughtful discussions with the members of the IFCC C-RSE.

11 References

- 1 ISO Guide 34:2009, General requirements for the competence of reference materials producers, International Organization for Standardization, Geneva, Switzerland.
- 2 ISO Guide 35:2006, Reference materials – General and statistical principles for certification, International Organization for Standardization, Geneva, Switzerland.
- 3 ISO/IEC 17025:2005, General requirements for the competence of testing and calibration laboratories, International Organization for Standardization, Geneva, Switzerland.
- 4 ISO/IEC Guide 98-3:2008, Guide to the expression of uncertainty in measurement, (GUM 1995), International Organization for Standardization, Geneva, Switzerland.
- 5 IUPAC-IUB Commission on Biochemical Nomenclature (CBN), Nomenclature of multiple forms of enzymes. Recommendations. *J. Biol. Chem.* 252 (1976) 5939-5941.
- 6 R.C. Haas, C. Korenfeld, Z. Zhang, B. Perryman, D. Roman, A.W. Strauss, Isolation and characterization of the gene cDNA encoding human mitochondrial CK. *J. Biol. Chem.* 264 (1989) 2890-2897.
- 7 N.A. Reiss, A.M. Kaye, Identification of the major component of the oestrogen-induced protein of rat uterus as the BB isoenzyme of creatine kinase. *J. Biol. Chem.* 256 (1981) 5741-5749.
- 8 L. Siekmann, R. Bonora, C.A. Burtis, F. Ceriotti, P. Clerc-Renaud, G. Féraud, C.A. Ferrero, J.-C. Forest, P.F.H. Franck, F.-J. Gella, W. Hoelzel, P.J. Jørgensen, T. Kanno, A. Kessner, IFCC primary reference procedures for the measurement of catalytic activity concentrations of enzymes at 37 °C - Part 1. The concept of reference procedures for the measurement of catalytic activity concentrations of enzymes, *Clin. Chem. Lab. Med.* 40 (2002) 631-634.
- 9 G. Schumann, R. Bonora, F. Ceriotti, P. Clerc-Renaud, C.A. Ferrero, G. Féraud, P.F.H. Franck, F.-J. Gella, W. Hoelzel, P.J. Jørgensen, T. Kanno, A. Kessner, R. Klauke, N. Kristiansen, J.-M. Lessinger, T.P.J. Linsinger, H. Misaki, M. Panteghini, J. Pauwels, H.G. Schimmel, A. Vialle, G. Weidemann, L. Siekmann, IFCC primary reference procedures for the measurement of catalytic activity concentrations of enzymes at 37 °C - Part 2. Reference procedure for the measurement of catalytic concentration of creatine kinase, *Clin. Chem. Lab. Med.* 40 (2002) 635-642.
- 10 G. Bonas, M. Zervou, T. Papaeoannou, M. Lees, "SoftCRM": A new software for the certification of reference materials, *Accred. Qual. Assur.*, 8 (2003) 101-107.
- 11 T.P.J. Linsinger, J. Pauwels, A.M.H. van der Veen, H. Schimmel, A. Lamberty, Homogeneity and stability of reference materials, *Accred. Qual. Assur.* 6 (2001) 20-25.
- 12 A. Lamberty, H. Schimmel, J. Pauwels, The study of the stability of reference materials by isochronous measurements, *Fres. J. Anal. Chem.* 360 (1998) 359-361.
- 13 T.P.J. Linsinger, J. Pauwels, A. Lamberty, H. Schimmel, A.M.H. van der Veen, L. Siekmann, Estimating the uncertainty of stability for matrix CRMs, *Fres. J. Anal. Chem.* 370 (2001) 183-188.
- 14 Comptes rendus de la 21e conférence générale des poids et mesures (1999), 2001, 334.

- 15 R. Bais, M. Panteghini, Chapter 15 Enzyme and Rate Analyses, In: Tietz textbook of Clinical Chemistry and Molecular Diagnostics, Editors: C. A. Burtis, E. R. Ashwood, D. E. Bruns, 5th Edition, 2011, 355-377.
- 16 R. Dybkær, Approved recommendation (1978) quantities and units in clinical chemistry, Clin. Chim. Acta 96 (1979) 157-183.
- 17 H. Vesper, H. Emons, M. Gnezda, C.P. Jain, W.G. Miller, R. Rej, G. Schumann, J. Tate, L. Thienpont, J.E. Vaks, Characterization and qualification of commutable reference materials for laboratory medicine; approved guideline, CLSI document C53-A, Clinical and Laboratory Standards Institute, Wayne, PA, USA, 2010.
- 18 T.P.J. Linsinger, ERM Application Note 1: Comparison of a measurement result with the certified value, www.erm-crm.org (last accessed on 20 May 2016).

Annexes

Annex A: Results of the homogeneity measurements

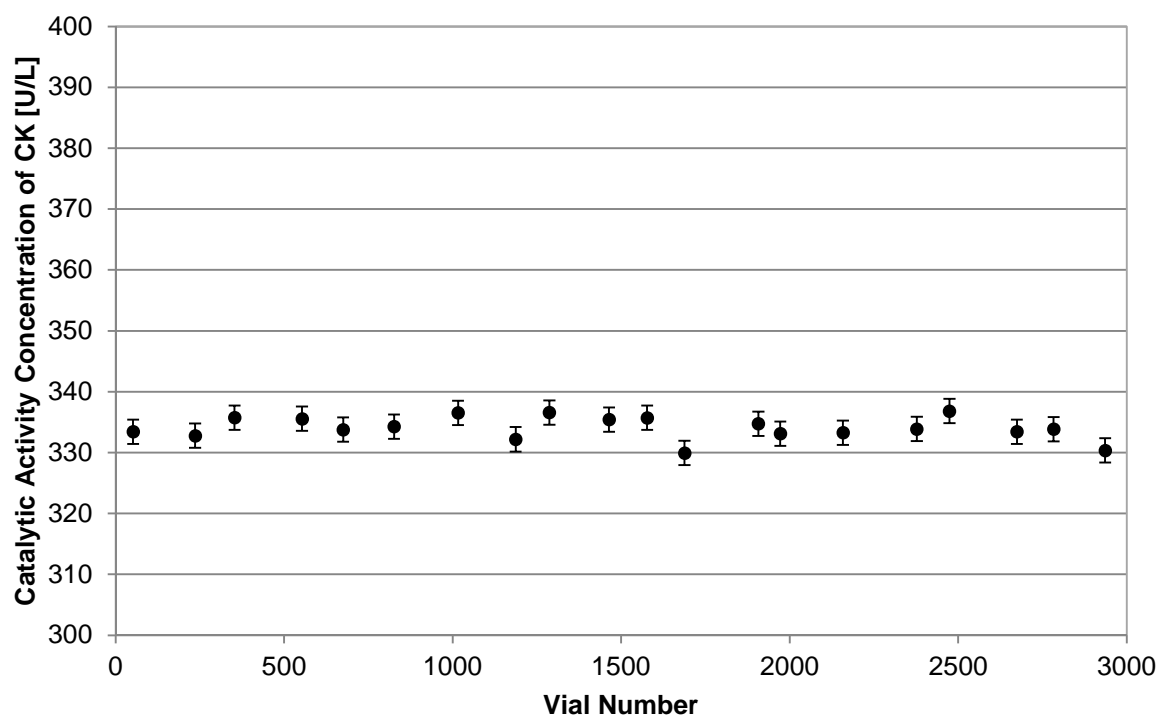


Figure A1: Homogeneity data of CK in ERM-AD455k/IFCC as measured with a UniCel DxC 800 Synchron Clinical System with CK reagent cartridges. Shown are the averages per vial number and their 95 % confidence interval based on the standard deviation as derived from a one-way ANOVA of all data.

Annex B: Results of the short-term stability measurements

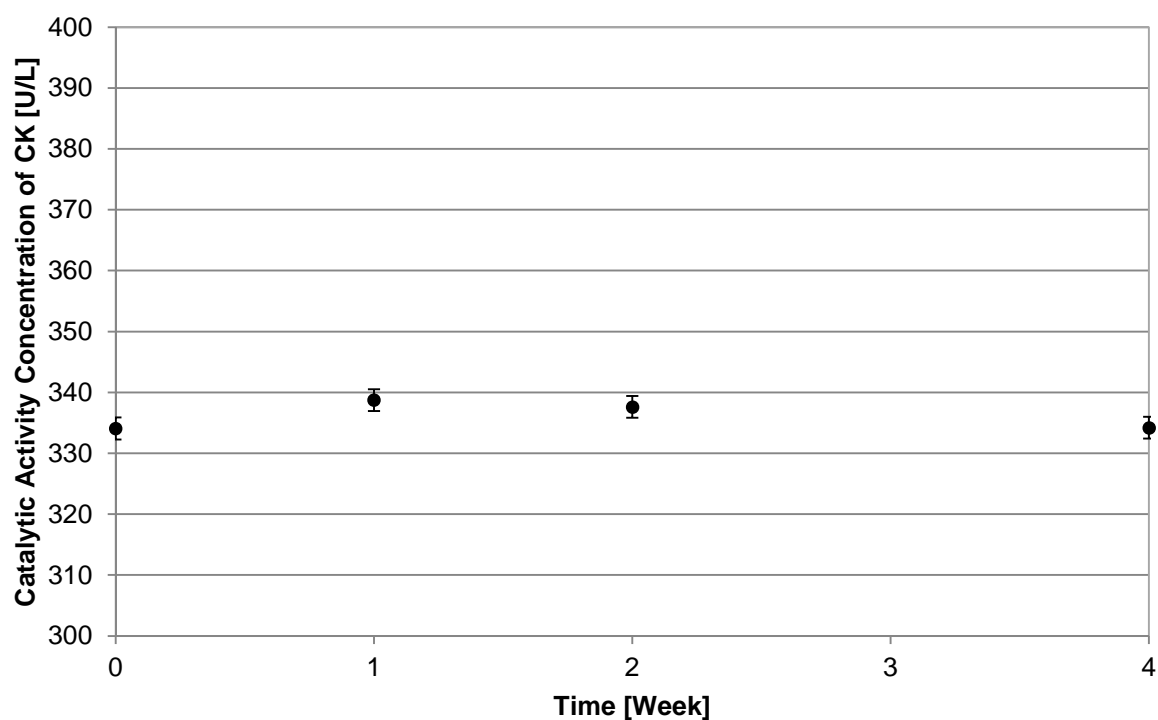


Figure B1: Short-term stability data of CK in ERM-AD455k/IFCC (stored at -20 °C) as measured with a UniCel DxC 800 Synchron Clinical System with CK reagent cartridges. Shown are the averages per time point and their 95 % confidence interval based on the standard deviation as derived from a one-way ANOVA of all data.

Annex C: Results of the long-term stability measurements

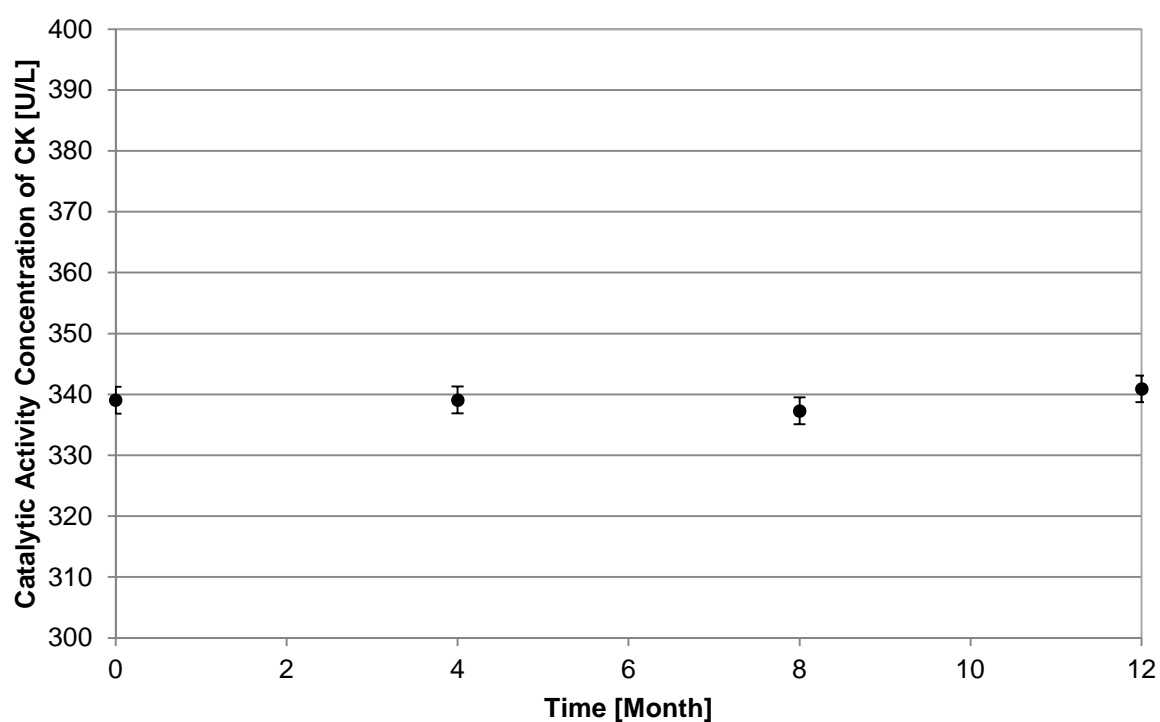


Figure C1: Long-term stability data of CK in ERM-AD455k/IFCC (stored at -20 °C) as measured with a UniCel DxC 800 Synchron Clinical System with CK reagent cartridges. Shown are the averages per time point and their 95 % confidence interval based on the standard deviation as derived from a one-way ANOVA of all data.

Annex D: Results of the characterisation measurements

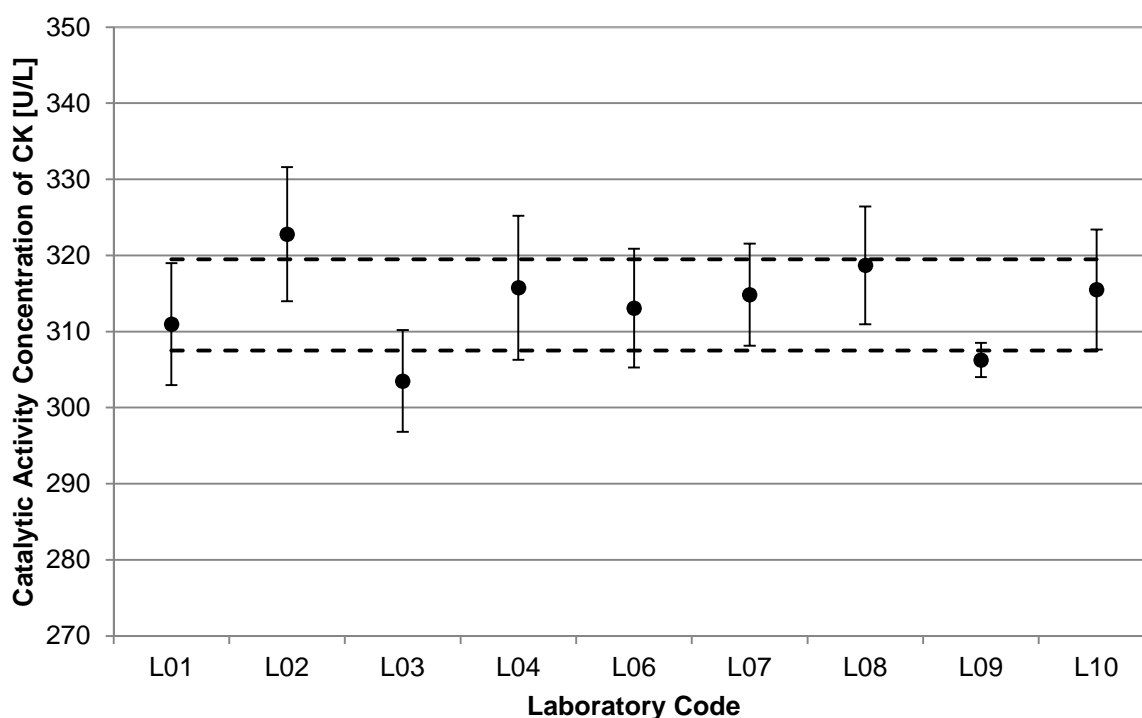


Figure D1: Graph showing average enzymatic activity concentrations of CK in ERM-AD455k/IFCC as measured with the IFCC primary reference measurement procedure at 37 °C with expanded uncertainties as stated by the laboratories and the 95 % certified interval (dotted lines).

Table D1: All accepted individual results, the mean value and the expanded uncertainty as stated by the laboratories for the catalytic activity concentration measurements of CK in ERM-AD455k/IFCC.

Laboratory code	Replicate 1 [U/L]	Replicate 2 [U/L]	Replicate 3 [U/L]	Replicate 4 [U/L]	Replicate 5 [U/L]	Replicate 6 [U/L]	Mean [U/L]	Expanded uncertainty [U/L]
L1	309.2	313.6	311.1	311.1	312.8	308.0	311.0	8.0
L2	322.74	321.77	324.40	321.41	324.86	321.63	322.8	8.8
L3	300	310	302	299	302	308	304	6.7
L4	318.3	313.8	317.2	316.6	315.2	313.4	315.8	9.47
L6	309.5	312.9	309.8	315.8	317.2	313.3	313.1	7.8
L7	317.76	313.60	313.72	316.01	313.93	313.97	314.83	6.7
L8	316.49	319.47	317.33	317.67	320.56	320.69	318.70	7.75
L9	307.1	306.3	304.5	307.4	306.9	305.4	306.3	2.237
L10	316	311	314	319	315	318	316	7.9

European Commission

EUR 28101 EN – Joint Research Centre – Directorate F – Health, Consumers and Reference Materials

Title: **The certification of the catalytic activity concentration of creatine kinase in
ERM®-AD455k/IFCC**

Author(s): Julia Kuhlmann, Brigitte Toussaint, Heinz Schimmel, Hendrik Emons

Luxembourg: Publications Office of the European Union

2016 – 26 pp. – 21.0 x 29.7 cm

EUR – Scientific and Technical Research series – ISSN 1831-9424

ISBN 978-92-79-61786-7

doi: 10.2787/310917

As the Commission's in-house science service, the Joint Research Centre's mission is to provide EU policies with independent, evidence-based scientific and technical support throughout the whole policy cycle.

Working in close cooperation with policy Directorates-General, the JRC addresses key societal challenges while stimulating innovation through developing new methods, tools and standards, and sharing its know-how with the Member States, the scientific community and international partners.

Key policy areas include: environment and climate change; energy and transport; agriculture and food security; health and consumer protection; information society and digital agenda; safety and security, including nuclear; all supported through a cross-cutting and multi-disciplinary approach.



doi: 10.2787/310917

ISBN: 978-92-79-61786-7