




Report on the 2007 Proficiency Test of the Community Reference Laboratory Network

Determination of Aflatoxins in a Peanut Product and a Test Solution

J. Stroka, A. Breidbach, I. Doncheva, M. Ambrosio

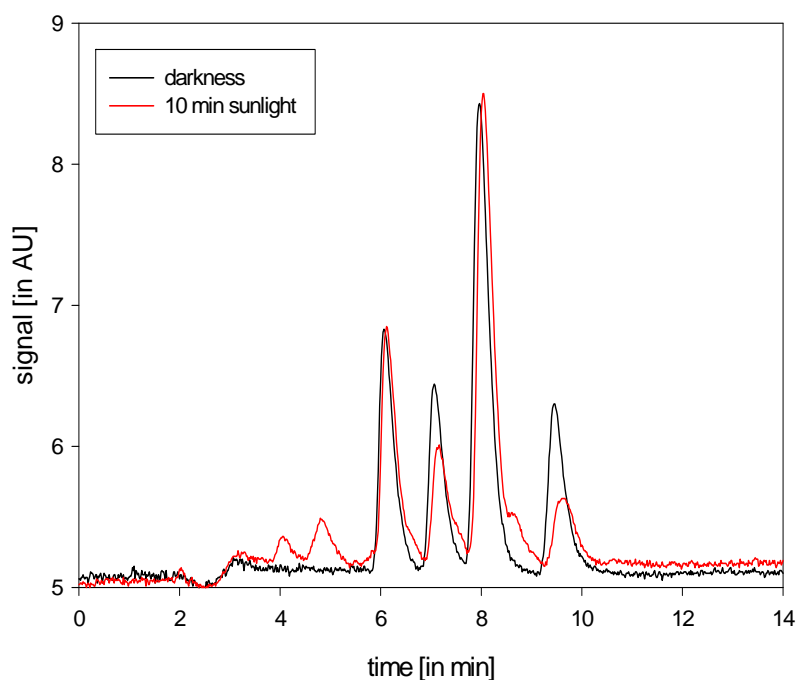

EUR 23261 EN - 2008

The mission of the IRMM is to promote a common and reliable European measurement system in support of EU policies.

European Commission
Joint Research Centre
Institute for Reference Materials and Measurements

Contact information

Address: Joerg Stroka, Retieseweg 111, B-2440 Geel
E-mail: Joerg.stroka@ec.europa.eu
Tel.: +32-14-571229
Fax: +32-14-571783

<http://irmm.jrc.ec.europa.eu/>
<http://www.jrc.ec.europa.eu/>

Legal Notice

Neither the European Commission nor any person acting on behalf of the Commission is responsible for the use which might be made of this publication.

***Europe Direct is a service to help you find answers
to your questions about the European Union***

Freephone number (*):

00 800 6 7 8 9 10 11

(*) Certain mobile telephone operators do not allow access to 00 800 numbers or these calls may be billed.

A great deal of additional information on the European Union is available on the Internet. It can be accessed through the Europa server <http://europa.eu/>

EUR 23261 EN
ISBN 978-92-79-08394-5
ISSN 1018-5593
DOI 10.2787/28654

Luxembourg: Office for Official Publications of the European Communities

© European Communities, 2008

Reproduction is authorised provided the source is acknowledged

Printed in Belgium

Table of Contents:

| Content | Page |
|------------------------------------|-------------|
| Summary..... | 4 |
| Introduction..... | 5 |
| Methodology..... | 5 |
| Results and Discussion..... | 6 |
| Conclusion..... | 17 |
| Annex..... | 18 |

Summary

A proficiency test was conducted with 31 European National Reference Laboratories (NRLs) for mycotoxins and one Laboratory from a candidate country. Test materials were a mixed aflatoxin (Af) solution in acetonitrile and two candidate Certified Reference Materials (CRM) - one "aflatoxin positive" and one "blank" material - that have not yet been released. Laboratories determined the aflatoxin content by reverse-phase high-performance liquid-chromatography (RP-HPLC) with either fluorescence or mass-selective detection against their own standard solutions as reference.

Applying the modified Horwitz equation according to Thompson¹ as a basis for the target standard deviation (22% in the case of this proficiency test), 26 out of 32 laboratories achieved z-scores of less than 2 and 17 laboratories reported values within the uncertainty range for both aflatoxin B1 and total aflatoxins in the candidate CRM after correction for recovery in both cases.

¹ M. Thompson (2000) *Analyst*, **125**, 385-386

Introduction

In 2006 the Institute for Reference Materials and Measurements (IRMM) in Geel was nominated as Community Reference Laboratory (CRL) for mycotoxins by the Directorate General for Health and Consumer Protection (DG SANCO). One of the main responsibilities of the CRL is to organise comparative testing to benchmark and harmonise the measurement capabilities of national reference laboratories (NRLs) working in the same field. Therefore, the CRL for mycotoxins together with the network of NRLs agreed to conduct the proficiency test in 2007 (PT2007) as follow up action to the PT 2006, this time on an aflatoxin solution in acetonitrile and a peanut material. This approach was chosen as it evaluates the first step, namely the calibration, as well as the analysis of the test material.

Methodology

Aflatoxin (Af) solutions were produced from the CRM calibrants that were used in the PT2006 and resulted in a mixed Af solution for which the assigned values are given in Table 1. In Table 2 the assigned values for the candidate CRM peanut material are given. The Af concentration values in the Af test solution was not known to the participants.

Table 1: Assigned values of the mixed aflatoxin test solution (the coverage factor $k=2$ corresponding to a level of confidence of about 95 %)

| Mixed Af solution | $\mu\text{g/L}$ | Uncertainty ($k=2$) |
|-------------------|-----------------|-----------------------|
| AfB1 | 20.85 | 0.61 |
| AfB2 | 5.95 | 0.13 |
| AfG1 | 5.90 | 0.21 |
| AfG2 | 5.95 | 0.12 |

Table 2: Assigned values of aflatoxins in peanut (the coverage factor $k=2$ corresponding to a level of confidence of about 95 %)

| Peanut CRM | $\mu\text{g/kg}$ | Uncertainty ($k=2$) |
|------------|------------------|-----------------------|
| AfB1 | 1.77 | 0.29 |
| AfB2 | 0.48 | 0.07 |
| AfG1 | 0.92 | 0.32 |
| AfG2 | 0.31 | 0.12 |

A full report on the production and certification of the CRM calibrants that were used to prepare the mixed aflatoxin calibrant solution is available from the IRMM. The report concerning the peanut CRMs will be available upon the final certification of these materials by the RM-Unit of the IRMM.

Each participant received an ampoule containing the Af test solution and two peanut materials ("Af positive" and "blank"). Participants were asked to measure the Af positive peanut material and the test solution for four aflatoxins, and to spike the blank peanut material with their own calibrant, reporting the spiking level and amount found to obtain recovery information. The instructions as sent to the participants are included in the annex.

Results and Discussion:

For each tested material the individual aflatoxin results are listed in **Tables 3 – 7** in the annex. As the repeatability for the three required measurements for each material and aflatoxin was in the lower %-range with an average of 5 %, all further calculations for the performance were made on the mean values calculated from the three measurements for each material/aflatoxin. In addition, no correlation was observed between the calculated repeatability values and the obtained z-scores. z-scores were calculated on the basis of the modified Horwitz equation according to Thompson. As a result, in all cases a target standard deviation of 22% was taken for $z = |1|$.

Deviation from the reference values

Figure 1 depicts the ranking of the results of the participating laboratories for aflatoxin B₁ in peanut prior recovery correction and prior calibrant correction. The order is by increasing laboratory mean value. The reference value and its uncertainty are depicted by a black and a blue line, respectively. The limit for a z-score of $z = |2|$ is indicated by red lines. For twelve of the laboratories the calculated laboratory mean values fell within that range of the assigned value. For twenty-six of the laboratories the calculated laboratory mean values fell within a z-score limit of 2. Seven laboratories reported values outside this limit.


Figure 2 depicts the ranking of the results of the participating laboratories for aflatoxin B₂ in peanut prior recovery correction and prior calibrant correction. The order is by increasing laboratory mean value. The reference value and its uncertainty are depicted by a black and a blue line, respectively. The limit for a z-score of $z = |2|$ is indicated by red lines for Figures 1 - 9. For thirteen of the laboratories the calculated laboratory mean values fell within that range of the assigned value. For twenty-four of the laboratories the calculated laboratory mean values fell within a z-score limit of 2.

Figure 3 depicts the ranking of the results of the participating laboratories for aflatoxin G₁ in peanut prior recovery correction and prior calibrant correction. The order is by increasing laboratory mean value. For twenty of the laboratories the calculated laboratory mean values fell within that range of the assigned value. For twenty-three of the laboratories the calculated laboratory mean values fell within a z-score limit of 2.

Figure 4 depicts the ranking of the results of the participating laboratories for aflatoxin G₂ in peanut prior recovery correction and prior calibrant correction. The order is by increasing laboratory mean value. For fourteen of the laboratories the calculated laboratory mean values fell within that range of the assigned value. For seventeen of the laboratories the calculated laboratory mean values fell within a z-score limit of 2.


Figure 5 depicts the ranking of the results of the participating laboratories for the sum of aflatoxins (total Af) in peanut prior recovery correction and prior calibrant correction. The order is by increasing laboratory mean value. For eleven of the laboratories the calculated laboratory mean values fell within that range of the assigned value. For twenty-six of the laboratories the calculated laboratory mean values fell within a z-score limit of 2.

Figure 1: Plot of Aflatoxin B1 in peanut prior recovery correction.

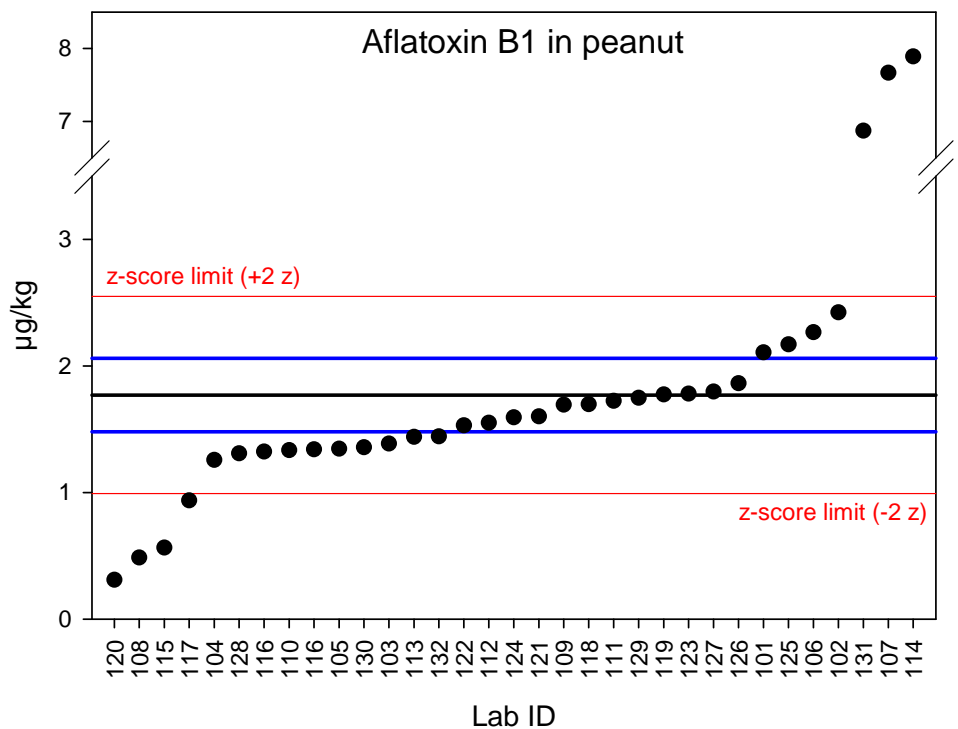


Plot of mean values from replicate measurements (n=3) for the determination of Aflatoxin B1 in a naturally contaminated peanut product. Blue lines reflect the uncertainty range of the reference value (black line), red lines the z-score limit of $z=|2|$.

Figure 2: Plot of Aflatoxin B2 in peanut prior recovery correction.

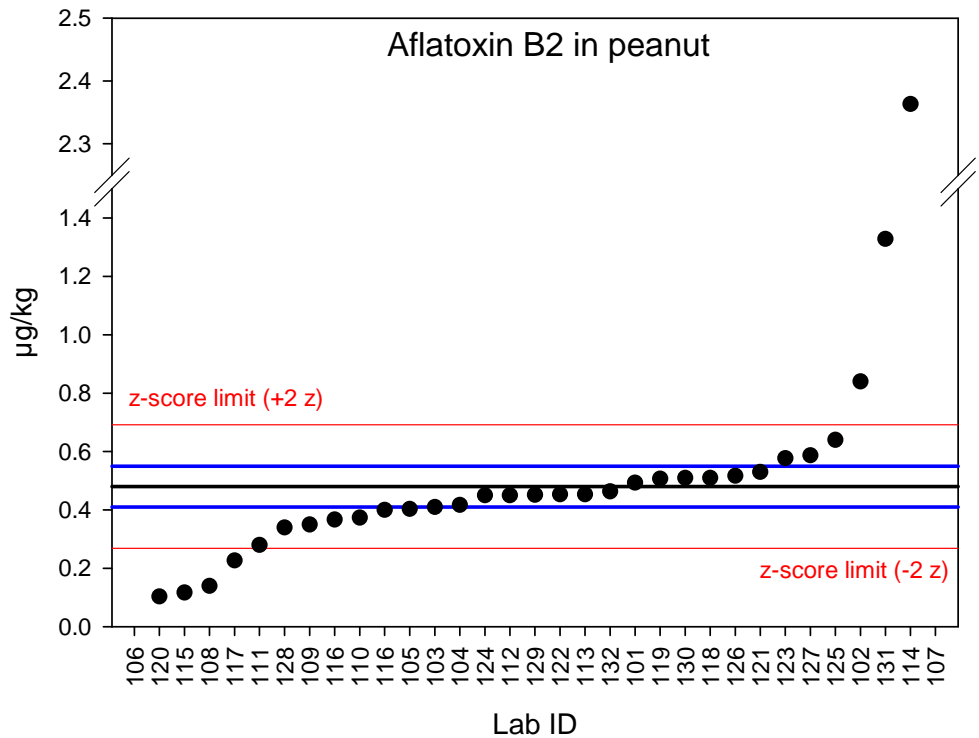


Plot of mean values from replicate measurements (n=3) for the determination of Aflatoxin B2 in a naturally contaminated peanut product. Blue lines reflect the uncertainty range of the reference value (black line), red lines the z-score limit of $z=|2|$.

Figure 3: Plot of Aflatoxin G1 in peanut prior recovery correction.

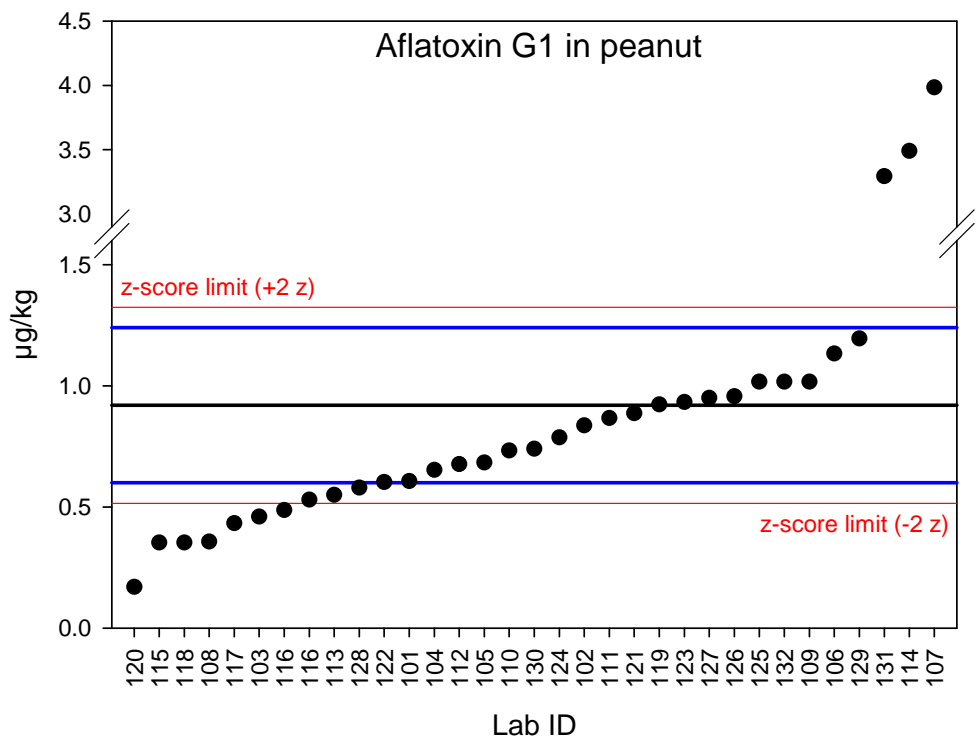


Plot of mean values from replicate measurements ($n=3$) for the determination of Aflatoxin G1 in a naturally contaminated peanut product. Blue lines reflect the uncertainty range of the reference value (black line), red lines the z-score limit of $z=|2|$.

Figure 4: Plot of Aflatoxin G2 in peanut prior recovery correction.

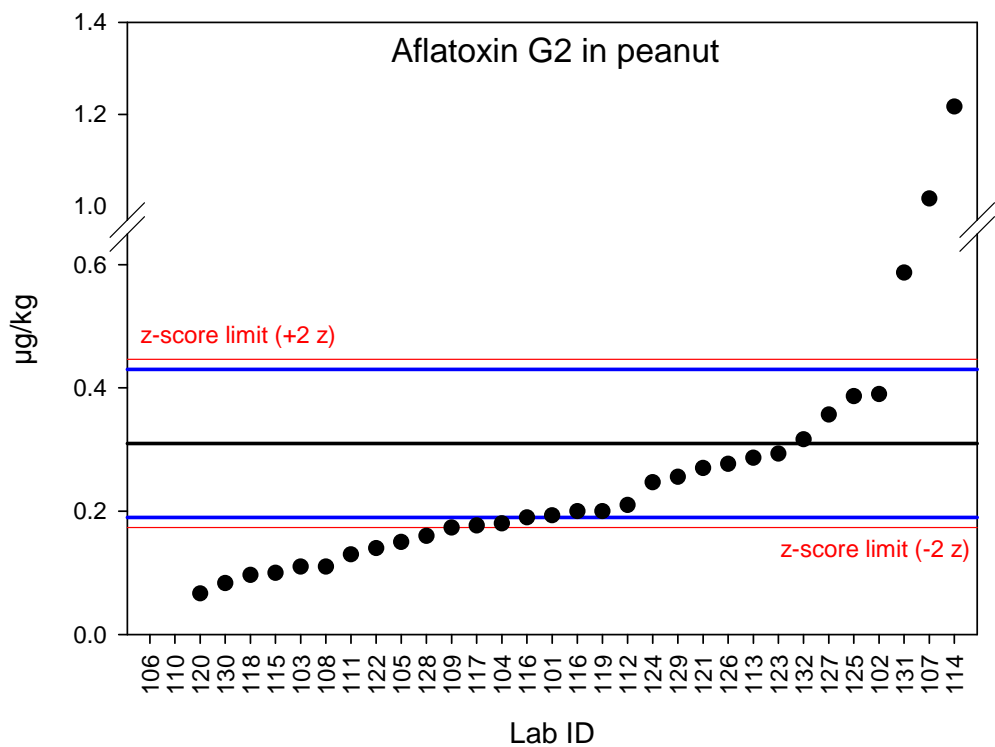

Plot of mean values from replicate measurements ($n=3$) for the determination of Aflatoxin G2 in a naturally contaminated peanut product. Blue lines reflect the uncertainty range of the reference value (black line), red lines the z-score limit of $z=|2|$.

Figure 5: Plot of total aflatoxins in peanut prior recovery correction.

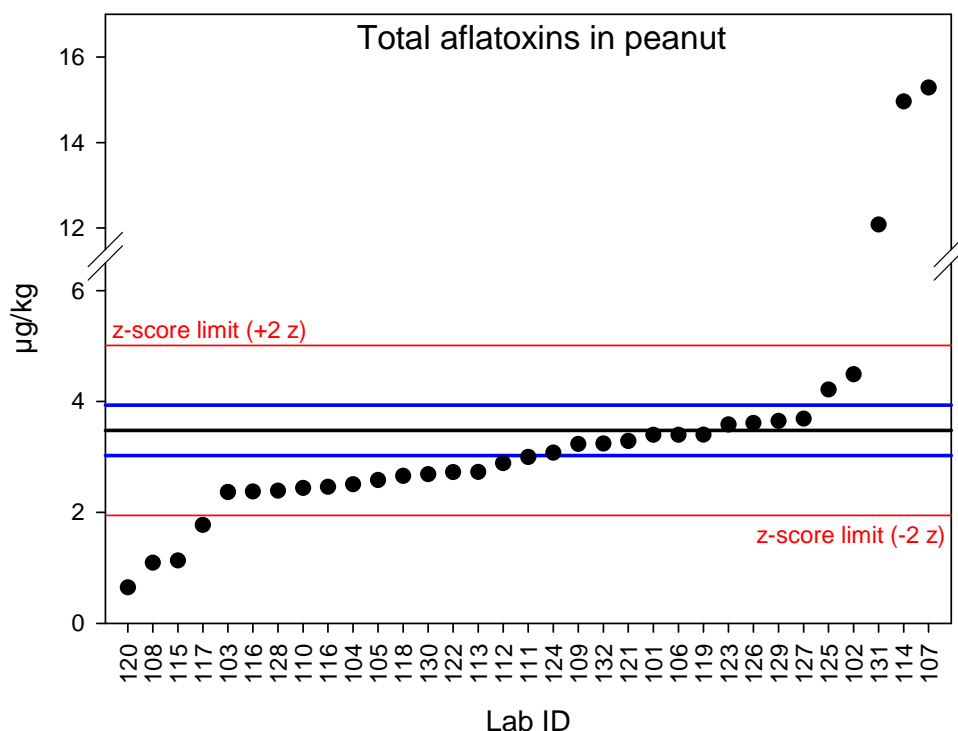


Plot of total aflatoxin values calculated from the sum of AfB1, AfB2, AfG1 and AfG2 in a naturally contaminated peanut product. Blue lines reflect the uncertainty range of the reference value (black line), red lines the z-score limit of $z=|2|$.

In addition to the analysis of a naturally contaminated sample ("aflatoxin positive") a blank peanut material was supplied. Participants were requested to indicate the level of aflatoxins spiked and analytically found, in order to calculate the recovery rate. Reported recoveries ranged from 53-118 % for AfB1, 56-142 % for AfB2, 46-106 % for AfG1 and 23-121 % for AfG2. The individually reported recovery figures were used to correct the results from the analysis of the "aflatoxin positive" material. For the determination of AfB1 and total Af results are shown in Figures 6 and 7. Table 3 summarises the effect by comparing the number of laboratories that reported values within the z-score limit and the uncertainty of the reference value.


Furthermore the results for AfB1 and total Af were corrected by a factor that was calculated from the measurement of the Af test solution. This factor was obtained by multiplying the result with the assigned value of the Af test solution, divided by the reported value (mean of 3 determinations). By this procedure values are normalised for calibrant effects. The results are shown in Figures 8 and 9.

Figure 6: Plot of Aflatoxin B1 results in peanut after recovery correction.

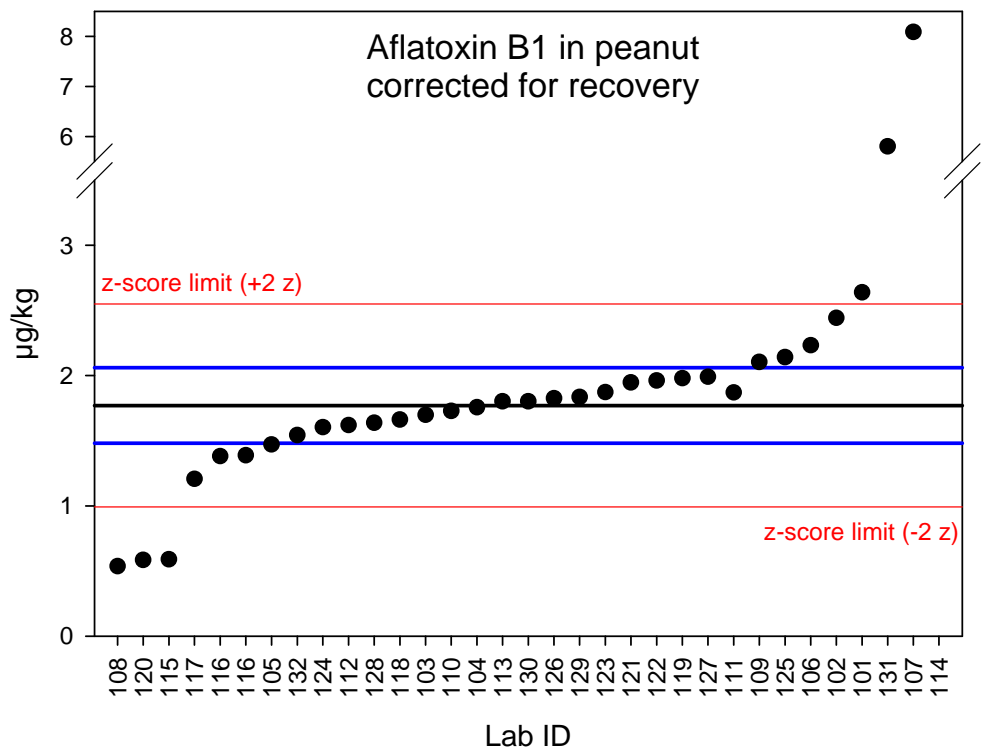


Plot of mean values from replicate measurements (n=3) for the determination of Aflatoxin B1 in a naturally contaminated peanut product after correction by recovery (mean of three replicate measurements). Blue lines reflect the uncertainty range of the reference value (black line), red lines the z-score limit of $z=|2|$.

Figure 7: Plot of total aflatoxins in peanut after recovery correction.

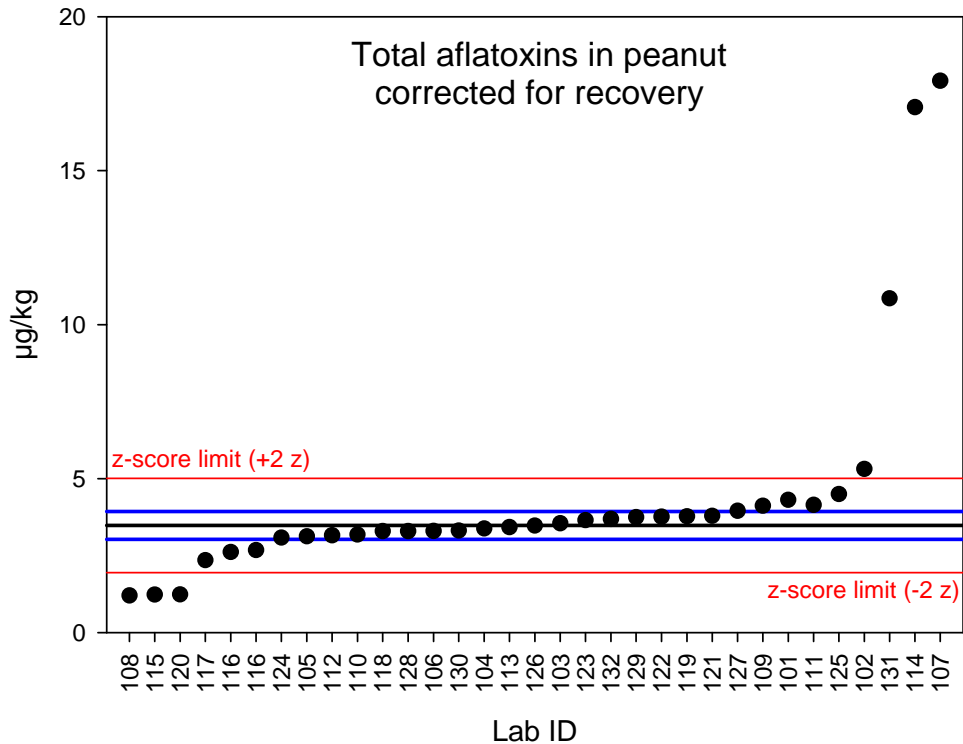


Plot of mean values from replicate measurements (n=3) for the sum of AfB1, AfB2, AfG1 and AfG2 in a naturally contaminated peanut product after correction by recovery (mean of three replicate measurements). Blue lines reflect the uncertainty range of the reference value (black line), red lines the z-score limit of $z=|2|$.

Figure 8: Plot of Aflatoxin B1 results in peanut after calibration bias correction.

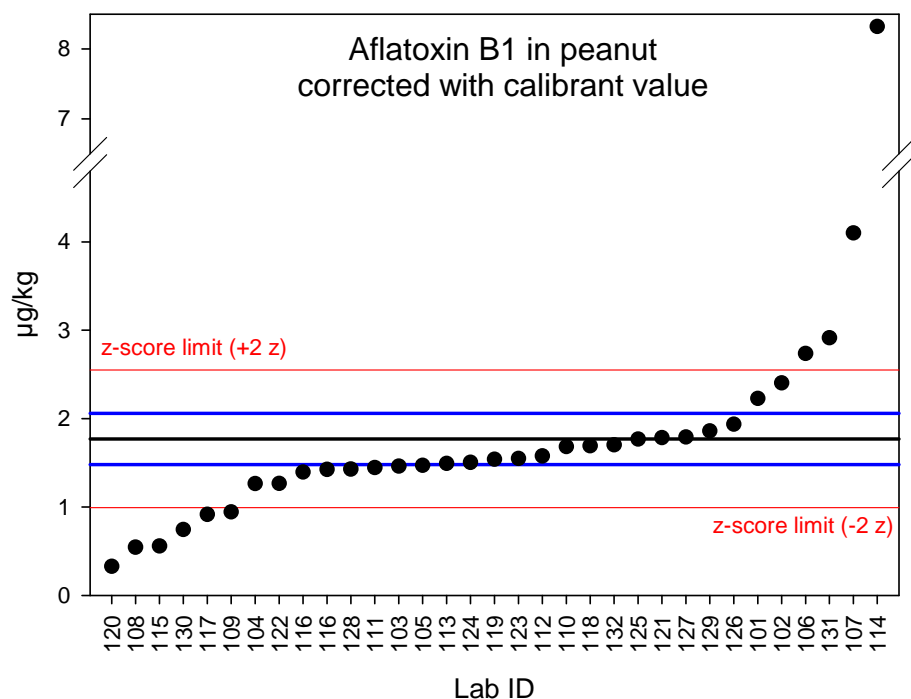

Plot of mean values from replicate measurements ($n=3$) for the determination of Aflatoxin B1 in a naturally contaminated peanut product after correction for the calibrant (mean of three replicate measurements). Blue lines reflect the uncertainty range of the reference value (black line), red lines the z-score limit of $z=|2|$.

Figure 9: Plot of total aflatoxins results in peanut after calibration bias correction.

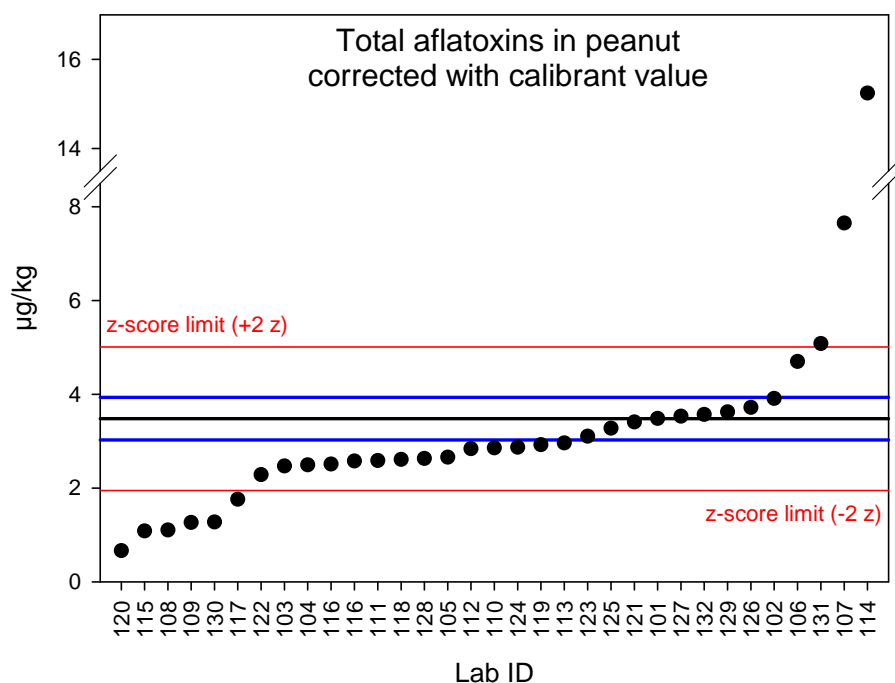

Plot of mean values from replicate measurements ($n=3$) for the sum of AfB1, AfB2, AfG1 and AfG2 in a naturally contaminated peanut product after correction for the calibrant (mean of three replicate measurements). Blue lines reflect the uncertainty range of the reference value (black line), red lines the z-score limit of $z=|2|$.

Table 3: Summary for values within z-score and certification limits

| n=33 | Number of results within z-score limit ($z \leq 2 $) | | | | | Number of results within certification limit | | | |
|-------|---|-----|---------|---------|---------|--|---------|---------|---------|
| | Af pos | Cal | Cal-cor | Rec-cor | Cal-rec | Af pos | Cal-cor | Rec-cor | Cal-rec |
| AfB1 | 26 | 29 | 23 | 26 | 25 | 12 | 13 | 17 | 17 |
| AfB2 | 24 | 29 | 23 | 25 | 25 | 13 | 15 | 17 | 15 |
| AfG1 | 23 | 29 | 21 | 26 | 25 | 20 | 17 | 21 | 23 |
| AfG2 | 17 | 28 | 17 | 20 | 22 | 14 | 14 | 20 | 21 |
| AfSUM | 26 | - | 24 | 26 | 26 | 11 | 9 | 18 | 16 |

Af pos = Af positive peanut prior any correction. Cal = mixed Af test solution. Cal-cor = Af pos after correction with the values obtained from the determination of the mixed Af test solution. Rec-cor = Af pos after recovery correction. Cal-rec = Af pos after calibrant and recovery correction.

Correction for recovery

As can be seen from **Figures 1, 5 - 7** and **Table 3** the correction for recovery has a clear effect on analytical results for the analysis of peanut material. This procedure aims to normalise the effects that occur during the analysis (extraction, clean-up, etc.). Thus, recovery correction improved slightly the number of acceptable z-scores and to a considerable extent the number of results within the uncertainty ranges of the reference values. This supports the finding that the correction for recovery improves the comparability of analytical results between laboratories². It should however be noted that the reported values which have been used for the correction were not obtained in a blinded experiment. Participants were asked to use their own calibrant, spike at a known level and report the found amount. From this the recovery was calculated.

Correction for calibrant

The intention to improve analytical results by normalising the effect that might be due to the quality of the calibrant was less effective than in the case of recovery or not effective at all. No improvement could be observed (see **Table 3**) as it was observed for recovery correction. An interesting observation is however that on an overall scale (all participants) the center-part of the ranked results in **Figure 8** is apparently more effected by the test solution correction than by recovery correction (**Figure 6**) for initial results of AfB1 (**Figure 1**). This is stressed by the fact that the "slope" of the plotted results is less steep in the middle region. This effect indicates that for a certain fraction of laboratories the test solution correction has a normalisation effect, too. Any further conclusion and evaluation would require a degree of data unscrambling that is beyond the scope of this proficiency test at this time.

Comparison of other observed trends within laboratories


Certain effects could be observed by plotting the reported AfB1 concentration values (as the major constituent in all test materials) found in relation to the assigned values of the three test materials (aflatoxin positive peanuts, peanut spike and Af test solution).

Ideally, all values reported by a laboratory should be located near the 100 % mark. Spike recovery for AfB1 as well as recovery of the Af test solution varied, with some exceptions, in a range between 70-110 %, while the AfB1 content reported in the aflatoxin positive material by the laboratories showed the sigmoidal distribution, which is typical for this type of PT exercise (**Figure 10**). Laboratories 120, 108, 115 and 117 had difficulty recovering the AfB1 from naturally incurred materials, although spike recovery was satisfactory; likewise laboratories 107 and 114 grossly overestimated the AfB1 content in the peanut material, while having satisfactory spike recovery. These findings suggest that for some experiments spiked or naturally incurred materials behaved differently during analysis, giving rise to an under or over estimation of the true content that can not be corrected using the results of the

² C. von Holst, J. Stroka, E. Anklam (2002), *Food Additives Contaminants*, **19**, 701-708

recovery experiment. Laboratories 130, 122 and 109 overestimated the AfB1 content in the AfB test solution, while reporting acceptable values for the natural and the spiked peanut materials. Incorrect dilution or handling of the AfB test solution might have caused these deviations.

Figure 10: Plot of AfB1 values in dependency of the type of analysis

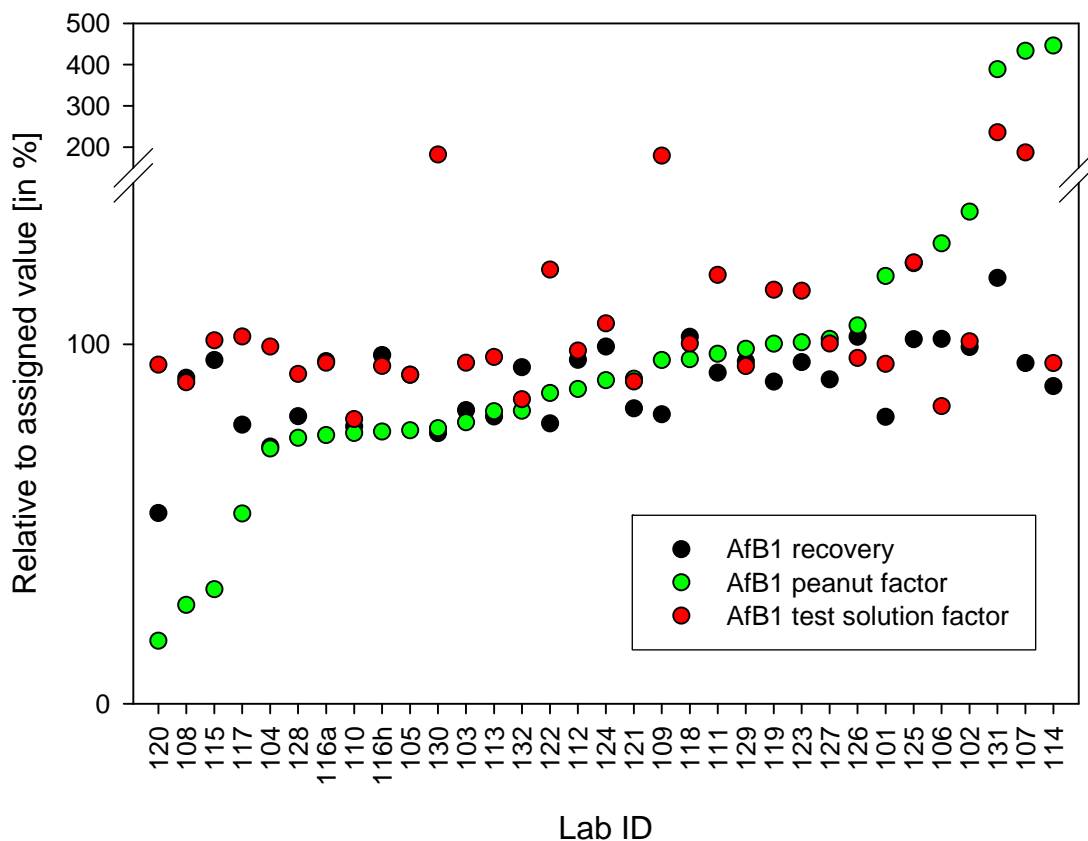

Values in the plot are ranked by the AfB1 value reported for the aflatoxin positive peanut material. Values were normalised to %-values of the respective assigned or fortified values.

Figure 11 shows that there is a general tendency to overestimate results for AfG2 in the Af test solution; this was not evident for the other aflatoxins. **Figure 12** illustrates that sample treatment (extraction, clean-up, etc.) increased the spread of Af values. However, the highest dispersion was found for AfG2 in this case with a tendency to lower values, which is in contrast to the findings for the Af test solution where generally higher values were found for AfG2. **Figure 13** shows that several laboratories reported considerable AfG2 losses in the recovery experiment, while for AfB1 the reported recoveries ranged from 80- 100 %. Nevertheless the low recovery values for AfG2 (and to a lesser extent also for AfG1) indicate a need for further investigation and corrective action, as such losses not necessarily occur at predictable rates and can thus lead to non-comparability and / or a misinterpretation of testing results. This is especially of relevance when AfG2 (and AfG1 – which is more likely) appear in higher concentrations in a material.

Figures 14 and 15 show the effectiveness of bias correction by using the Af testing solution. Based on the assumption that the recovery rate was properly established such a plot visualises effects related to the *in-house* calibrant (of *known* concentration) used by the laboratories. In cases where the values for the test solution results differ to a larger extent from the recovery corrected results of the Af, an *in-house* calibrant effect (=bias) is very likely. In those cases where the recovery corrected values for the peanut material were different from values of the Af test solution (values for the test solution being near 100%), a doubtful spike recovery estimate or problems in the course of the analytical procedure are likely.

These scenarios must be judged on the basis that the levels of aflatoxins G1 and G2 in the peanut material were rather low compared to levels for AfB1 and AfB2. This fact makes an exact evaluation/unscrambling of effects difficult for the G aflatoxins. Furthermore legislative limits exist currently only for AfB1 and total aflatoxins and thus these should be the final parameters focussed on. This however shall not mean that analytical performance for the less abundant aflatoxins can be neglected. An evaluation of effects from analytical procedures based on the questionnaire has not been performed, but all answers have been added to the annex to allow further unscrambling of data by the participants if necessary in particular cases.

Figure 11: Plot of the normalised aflatoxin values reported for the calibrant

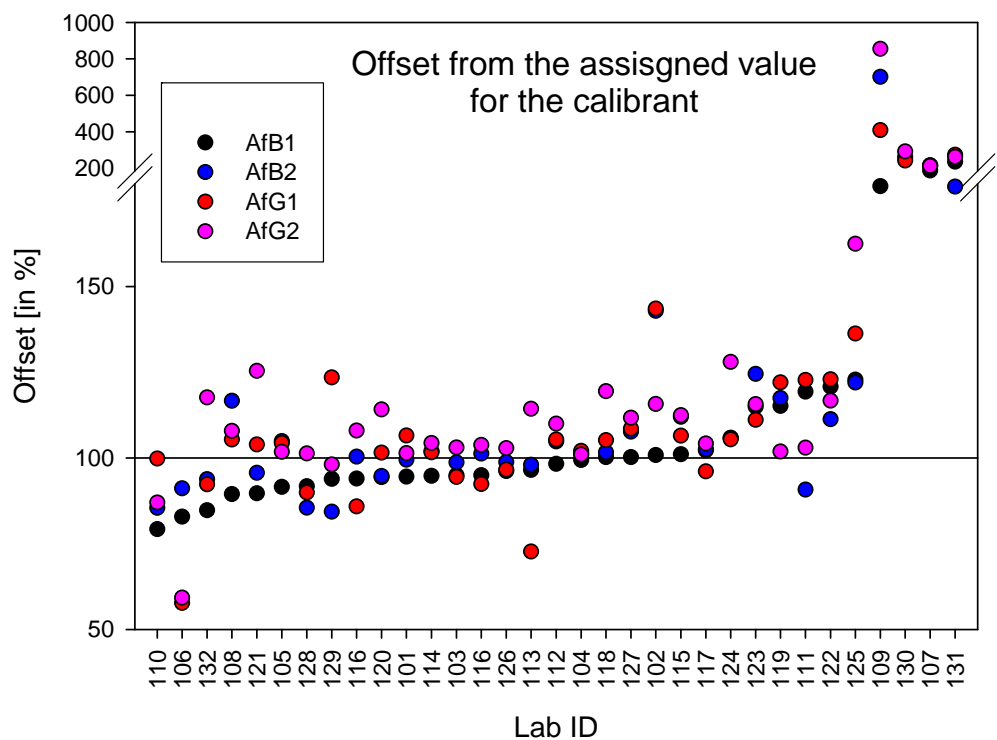


Figure 12: Plot of the normalised aflatoxin values reported for the peanut material

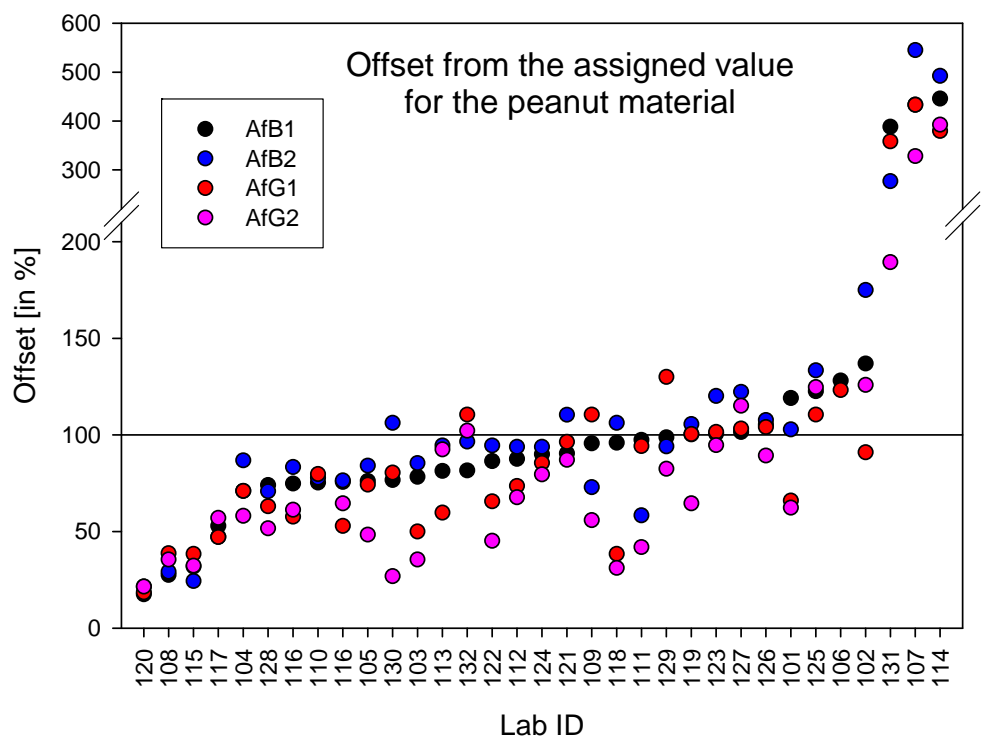


Figure 13: Plot of the aflatoxin recovery values reported for the spike

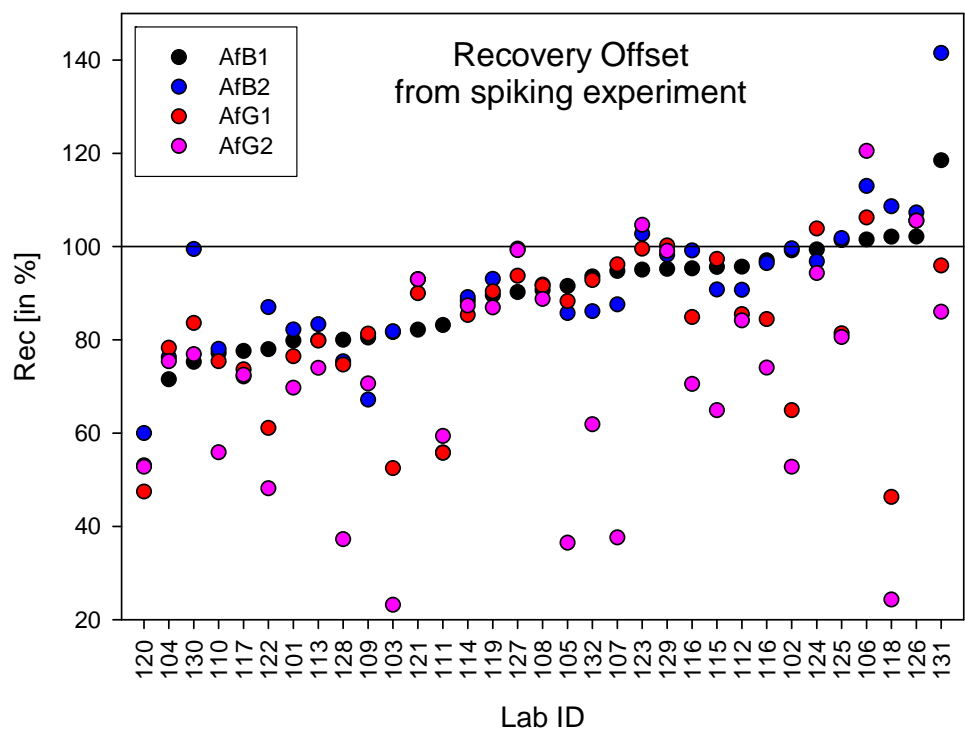


Figure 14: Plot of the normalised calibrant values for AfB1 and AfB2 against the recovery corrected result for the peanut positive result


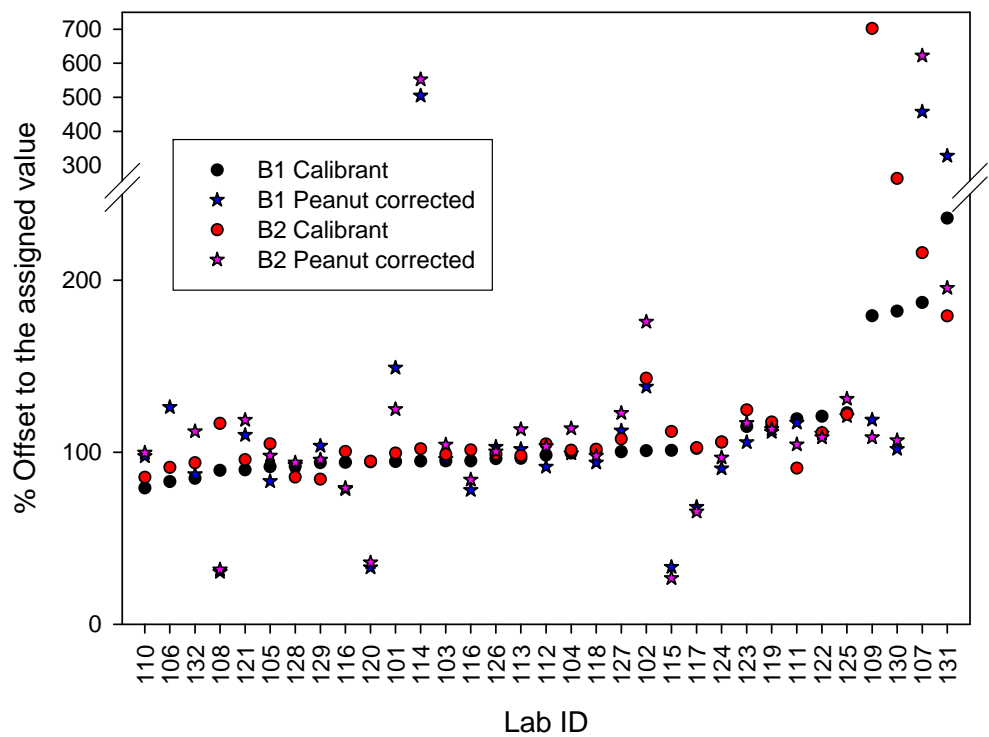


Figure 15: Plot of the normalised calibrant values for AfG1 and AfG2 against the recovery corrected result for the peanut positive result

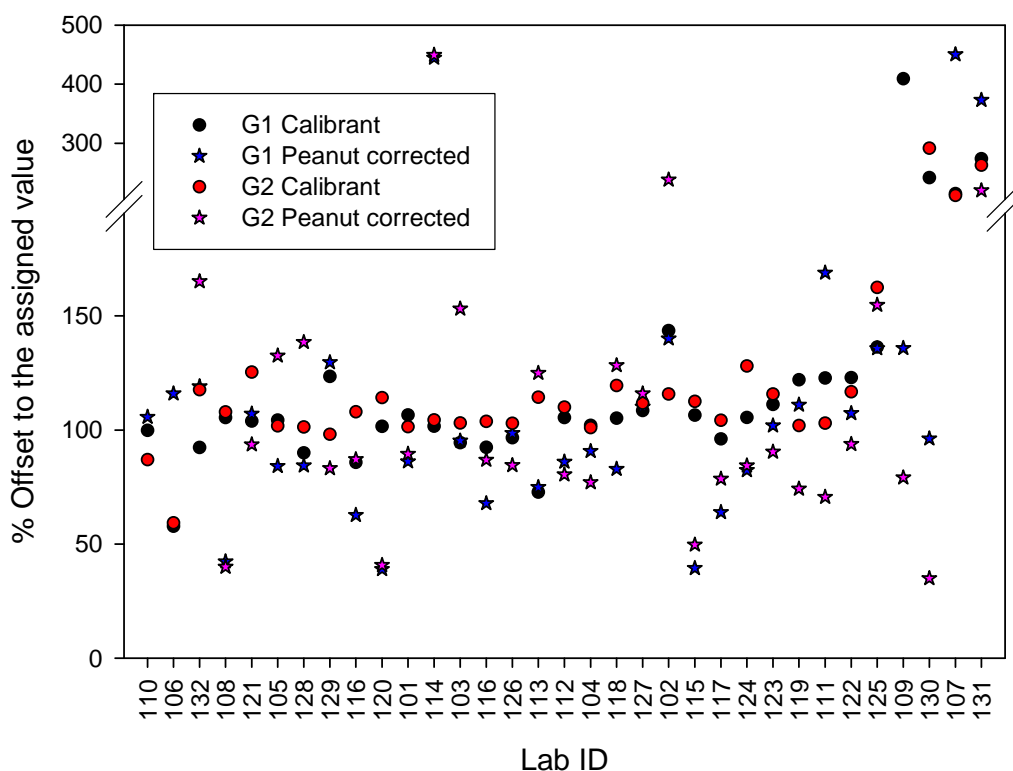

Conclusion

With respect to the NRL proficiency test 2006 the situation has improved. The majority of laboratories reported values within z-scores of $|2|$ taking the modified Horwitz equation as basis for the target standard deviation. For measurements of Af in a test solution 78% of the labs were able to report acceptable values for aflatoxin B1, while for aflatoxin G2 this was only the case for 59% of the laboratories. In general, the agreement of results among laboratories as well as in relation to the assigned values - improved after recovery correction. A further correction of results for AfB1 and total aflatoxins with a factor to take into account the calibration bias had no significant effect on the overall comparability of such corrected results, while a slight improvement for generally good performing labs could be observed. This indicates that a calibration bias correction can have a positive effect, however not in all cases.

Annex

Table 3: Individual results [in µg/kg] for aflatoxin B1 and B2 in peanut:

| Laboratory code | AfB1 | AfB1 | AfB1 | AfB2 | AfB2 | AfB2 |
|-----------------|-------|-------|-------|-------|-------|-------|
| 101 | 2.07 | 2.11 | 2.14 | 0.49 | 0.49 | 0.5 |
| 102 | 2.65 | 2.21 | 2.41 | 0.71 | 0.86 | 0.95 |
| 103 | 1.39 | 1.39 | 1.38 | 0.41 | 0.41 | 0.41 |
| 104 | 1.28 | 1.28 | 1.21 | 0.38 | 0.37 | 0.5 |
| 105 | 1.39 | 1.28 | 1.37 | 0.43 | 0.39 | 0.39 |
| 106 | 2.1 | 2.4 | 2.3 | <1 | <1 | <1 |
| 107 | 8.15 | 7.15 | 7.7 | 2.35 | 2.55 | 2.95 |
| 108 | 0.5 | 0.48 | 0.48 | 0.14 | 0.14 | 0.14 |
| 109 | 1.94 | 1.8 | 1.34 | 0.39 | 0.39 | 0.27 |
| 110 | 1.33 | 1.34 | 1.33 | 0.38 | 0.37 | 0.37 |
| 111 | 1.74 | 1.67 | 1.76 | 0.29 | 0.27 | 0.28 |
| 112 | 1.56 | 1.53 | 1.56 | 0.45 | 0.45 | 0.45 |
| 113 | 1.44 | 1.43 | 1.45 | 0.45 | 0.45 | 0.46 |
| 114 | 7.98 | 7.67 | 8.02 | 2.33 | 2.29 | 2.47 |
| 115 | 0.58 | 0.56 | 0.55 | 0.12 | 0.12 | 0.11 |
| 116 | 1.33 | 1.33 | 1.36 | 0.4 | 0.4 | 0.4 |
| 116 | 1.31 | 1.33 | 1.33 | 0.36 | 0.36 | 0.38 |
| 117 | 0.93 | 0.93 | 0.95 | 0.23 | 0.22 | 0.23 |
| 118 | 1.83 | 1.64 | 1.62 | 0.56 | 0.49 | 0.48 |
| 119 | 1.7 | 1.83 | 1.79 | 0.49 | 0.52 | 0.51 |
| 120 | 0.28 | 0.31 | 0.34 | 0.08 | 0.12 | 0.11 |
| 121 | 1.6 | 1.6 | 1.6 | 0.53 | 0.53 | 0.53 |
| 122 | 1.49 | 1.59 | 1.51 | 0.45 | 0.45 | 0.46 |
| 123 | 1.67 | 1.8 | 1.87 | 0.53 | 0.59 | 0.61 |
| 124 | 1.56 | 1.55 | 1.67 | 0.46 | 0.45 | 0.44 |
| 125 | 2.27 | 2.08 | 2.16 | 0.64 | 0.63 | 0.65 |
| 126 | 1.82 | 1.86 | 1.91 | 0.49 | 0.53 | 0.53 |
| 127 | 1.79 | 1.8 | 1.8 | 0.59 | 0.58 | 0.59 |
| 128 | 1.3 | 1.28 | 1.35 | 0.34 | 0.34 | 0.34 |
| 129 | 1.722 | 1.774 | 1.746 | 0.452 | 0.456 | 0.446 |
| 130 | 1.46 | 1.3 | 1.31 | 0.57 | 0.48 | 0.48 |
| 131 | 6.921 | 6.893 | 6.806 | 1.331 | 1.337 | 1.315 |
| 132 | 1.41 | 1.47 | 1.45 | 0.46 | 0.47 | 0.46 |

Table 4: Individual results [in µg/kg] for aflatoxin G1 and G2 in peanut:

| Laboratory code | AfG1 | AfG1 | AfG1 | AfG2 | AfG2 | AfG2 |
|-----------------|-------|-------|------|-------|-------|-------|
| 101 | 0.63 | 0.61 | 0.58 | 0.19 | 0.2 | 0.19 |
| 102 | 0.8 | 0.82 | 0.89 | 0.39 | 0.39 | 0.39 |
| 103 | 0.46 | 0.46 | 0.46 | 0.11 | 0.11 | 0.11 |
| 104 | 0.74 | 0.6 | 0.62 | 0.19 | 0.18 | 0.17 |
| 105 | 0.7 | 0.66 | 0.69 | 0.16 | 0.16 | 0.13 |
| 106 | 1.2 | 1 | 1.2 | <1 | <1 | <1 |
| 107 | 4.1 | 3.7 | 4.15 | 1 | 1.25 | 0.8 |
| 108 | 0.35 | 0.35 | 0.37 | 0.11 | 0.11 | 0.11 |
| 109 | 1.11 | 1.01 | 0.93 | 0.19 | 0.18 | 0.15 |
| 110 | 0.75 | 0.74 | 0.71 | <0.21 | <0.21 | <0.21 |
| 111 | 0.85 | 0.87 | 0.88 | 0.14 | 0.12 | 0.13 |
| 112 | 0.68 | 0.68 | 0.67 | 0.2 | 0.21 | 0.22 |
| 113 | 0.55 | 0.55 | 0.55 | 0.29 | 0.28 | 0.29 |
| 114 | 3.71 | 3.42 | 3.34 | 1.39 | 1.08 | 1.18 |
| 115 | 0.35 | 0.35 | 0.36 | 0.1 | 0.1 | 0.1 |
| 116 | 0.53 | 0.53 | 0.53 | 0.19 | 0.19 | 0.19 |
| 116 | 0.49 | 0.48 | 0.49 | 0.2 | 0.19 | 0.21 |
| 117 | 0.43 | 0.43 | 0.44 | 0.18 | 0.17 | 0.18 |
| 118 | 0.31 | 0.39 | 0.36 | 0.09 | 0.09 | 0.11 |
| 119 | 0.88 | 0.94 | 0.95 | 0.2 | 0.22 | 0.18 |
| 120 | 0.17 | 0.18 | 0.16 | 0.06 | 0.08 | 0.06 |
| 121 | 0.93 | 0.93 | 0.8 | 0.27 | 0.27 | 0.27 |
| 122 | 0.6 | 0.58 | 0.63 | 0.15 | 0.13 | 0.14 |
| 123 | 0.86 | 0.96 | 0.98 | 0.27 | 0.3 | 0.31 |
| 124 | 0.79 | 0.79 | 0.78 | 0.24 | 0.24 | 0.26 |
| 125 | 1.03 | 0.96 | 1.06 | 0.41 | 0.35 | 0.4 |
| 126 | 0.87 | 1.02 | 0.98 | 0.23 | 0.31 | 0.29 |
| 127 | 0.95 | 0.95 | 0.95 | 0.36 | 0.36 | 0.35 |
| 128 | 0.59 | 0.57 | 0.58 | 0.16 | 0.16 | 0.16 |
| 129 | 1.216 | 1.181 | 1.19 | 0.234 | 0.259 | 0.274 |
| 130 | 0.88 | 0.67 | 0.67 | 0.08 | 0.1 | 0.07 |
| 131 | 3.294 | 3.306 | 3.28 | 0.628 | 0.621 | 0.512 |
| 132 | 1.03 | 0.99 | 1.03 | 0.39 | 0.24 | 0.32 |

Table 5: Individual results [in µg/kg] for aflatoxin B1 and B2 in the aflatoxin test solution:

| Laboratory code | AfB1 | AfB1 | AfB1 | AfB2 | AfB2 | AfB2 |
|-----------------|--------|--------|--------|-------|--------|--------|
| 101 | 19.4 | 19.78 | 19.8 | 5.86 | 5.89 | 5.95 |
| 102 | 21.48 | 21.07 | 20.33 | 8.87 | 8.78 | 7.77 |
| 103 | 20.06 | 19.72 | 19.42 | 6.04 | 5.82 | 5.7 |
| 104 | 21 | 21 | 20 | 6 | 6 | 6 |
| 105 | 18.98 | 19.01 | 19.1 | 6.28 | 6.28 | 6.1 |
| 106 | 16.43 | 17.188 | 18.044 | 5.425 | 5.727 | 5.059 |
| 107 | 39.5 | 38.54 | 38.59 | 12.99 | 12.62 | 12.8 |
| 108 | 18.514 | 18.631 | 18.627 | 6.893 | 6.928 | 6.939 |
| 109 | 39.85 | 35.18 | 36.79 | 34.94 | 40.07 | 49.87 |
| 110 | 16.3 | 16.9 | 16.2 | 5 | 5.1 | 5.1 |
| 111 | 24.89 | 24.85 | 24.69 | 5.43 | 5.384 | 5.331 |
| 112 | 20.381 | 20.413 | 20.503 | 6.255 | 6.198 | 6.187 |
| 113 | 19.94 | 19.98 | 20.25 | 5.81 | 5.78 | 5.84 |
| 114 | 19.919 | 19.786 | 19.415 | 6.066 | 6.125 | 5.955 |
| 115 | 21.11 | 20.86 | 21.08 | 6.64 | 6.63 | 6.65 |
| 116 | 19.23 | 19.69 | 19.71 | 5.86 | 6.12 | 6.03 |
| 116 | 20.12 | 19.64 | 19.44 | 6.04 | 5.97 | 5.85 |
| 117 | 21.689 | 21.303 | 20.732 | 6.167 | 6.058 | 6.023 |
| 118 | 20.82 | 20.97 | 20.73 | 6 | 6.12 | 5.96 |
| 119 | 24.21 | 23.97 | 23.67 | 7.08 | 6.85 | 6.97 |
| 120 | 19.799 | 19.235 | 19.823 | 5.667 | 5.481 | 5.703 |
| 121 | 18.73 | 18.6 | 18.6 | 5.67 | 5.67 | 5.67 |
| 122 | 25.32 | 25.02 | 25 | 6.58 | 6.52 | 6.7 |
| 123 | 23.685 | 23.999 | 24.002 | 7.388 | 7.394 | 7.37 |
| 124 | 22.001 | 22.003 | 22.016 | 6.276 | 6.277 | 6.278 |
| 125 | 24.86 | 25.18 | 26.55 | 6.95 | 7.1 | 7.65 |
| 126 | 19.286 | 19.951 | 20.77 | 5.708 | 5.859 | 6.011 |
| 127 | 20.835 | 20.665 | 21.021 | 6.363 | 6.31 | 6.475 |
| 128 | 18.27 | 18.68 | 20.26 | 4.98 | 5 | 5.23 |
| 129 | 19.41 | 19.64 | 19.52 | 5 | 5.06 | 4.94 |
| 130 | 37.93 | 37.95 | 37.7 | 15.71 | 15.48 | 15.28 |
| 131 | 48.262 | 49.117 | 49.82 | 10.09 | 10.807 | 10.972 |
| 132 | 17.669 | 17.56 | 17.628 | 5.633 | 5.503 | 5.551 |

Table 6: Individual results [in µg/kg] for aflatoxin G1 and G2 in the aflatoxin test solution:

| Laboratory code | AfG1 | AfG1 | AfG1 | AfG2 | AfG2 | AfG2 |
|-----------------|-------|--------|--------|--------|--------|-------|
| 101 | 6.14 | 6.41 | 6.24 | 5.95 | 6.07 | 6.05 |
| 102 | 9.03 | 8.09 | 8.2 | 7.63 | 7.08 | 5.91 |
| 103 | 5.72 | 5.54 | 5.4 | 6.16 | 6.12 | 6.08 |
| 104 | 6 | 6 | 6 | 6 | 6 | 6 |
| 105 | 6.18 | 6.17 | 6.05 | 6.16 | 6.14 | 5.84 |
| 106 | 3.131 | 3.063 | 3.989 | 3.48 | 3.711 | 3.373 |
| 107 | 12.75 | 12.65 | 12.42 | 12.77 | 12.64 | 12.27 |
| 108 | 6.171 | 6.211 | 6.211 | 6.383 | 6.418 | 6.425 |
| 109 | 29.64 | 22.62 | 19.9 | 48.27 | 51.06 | 52.98 |
| 110 | 5.8 | 6 | 5.8 | 5.6 | 4.8 | 5.1 |
| 111 | 7.288 | 7.137 | 7.225 | 6.148 | 6.11 | 6.093 |
| 112 | 6.433 | 5.919 | 6.24 | 6.832 | 6.301 | 6.465 |
| 113 | 4.2 | 4.2 | 4.43 | 6.79 | 6.7 | 6.88 |
| 114 | 6.015 | 6.07 | 5.84 | 6.241 | 6.302 | 6.049 |
| 115 | 6.26 | 6.3 | 6.22 | 6.71 | 6.74 | 6.59 |
| 116 | 5.35 | 5.62 | 5.33 | 6.13 | 6.28 | 6.08 |
| 116 | 5.11 | 5.13 | 4.9 | 6.51 | 6.47 | 6.26 |
| 117 | 5.647 | 5.626 | 5.673 | 6.213 | 6.171 | 6.182 |
| 118 | 6.19 | 6.29 | 6.07 | 7.15 | 7.13 | 7.01 |
| 119 | 7.36 | 7.15 | 7.01 | 6.25 | 5.83 | 6.08 |
| 120 | 6.041 | 5.889 | 5.978 | 6.812 | 6.545 | 6.981 |
| 121 | 6.13 | 6.13 | 6.07 | 7.47 | 7.47 | 7.4 |
| 122 | 7.16 | 7.22 | 7.3 | 6.84 | 6.9 | 7.06 |
| 123 | 6.562 | 6.544 | 6.497 | 6.87 | 6.903 | 6.843 |
| 124 | 6.203 | 6.193 | 6.203 | 7.603 | 7.602 | 7.604 |
| 125 | 8.03 | 7.86 | 8.15 | 9.54 | 9.46 | 9.95 |
| 126 | 5.256 | 5.779 | 6.007 | 5.851 | 6.246 | 6.23 |
| 127 | 6.56 | 6.355 | 6.228 | 6.608 | 6.697 | 6.601 |
| 128 | 5.16 | 5.25 | 5.45 | 5.98 | 6.01 | 6.06 |
| 129 | 7.07 | 7.51 | 7.2 | 5.82 | 5.86 | 5.8 |
| 130 | 14.39 | 14.54 | 13.71 | 17.57 | 17.24 | 17.14 |
| 131 | 16.09 | 16.242 | 15.894 | 15.672 | 15.575 | 15.58 |
| 132 | 5.496 | 5.362 | 5.42 | 6.907 | 6.978 | 7.074 |

Table 7: Mean spike recoveries calculated for aflatoxin B1, B2, G1 and G2 [in %]:

| Laboratory code | AfB1 | AfB1 | AfG1 | AfG2 |
|-----------------|------|------|------|------|
| 101 | 80 | 82 | 77 | 70 |
| 102 | 99 | 100 | 65 | 53 |
| 103 | 82 | 82 | 53 | 23 |
| 104 | 72 | 76 | 78 | 75 |
| 105 | 92 | 86 | 88 | 37 |
| 106 | 102 | 113 | 106 | 121 |
| 107 | 95 | 88 | 96 | 38 |
| 108 | 91 | 92 | 92 | 89 |
| 109 | 81 | 67 | 81 | 71 |
| 110 | 77 | 78 | 75 | 56 |
| 111 | 92 | 56 | 56 | 59 |
| 112 | 96 | 91 | 86 | 84 |
| 113 | 80 | 83 | 80 | 74 |
| 114 | 88 | 89 | 85 | 87 |
| 115 | 96 | 91 | 97 | 65 |
| 116 | 97 | 99 | 85 | 71 |
| 116 | 95 | 96 | 84 | 74 |
| 117 | 78 | 72 | 74 | 72 |
| 118 | 102 | 109 | 46 | 24 |
| 119 | 90 | 93 | 90 | 87 |
| 120 | 53 | 60 | 47 | 53 |
| 121 | 82 | 93 | 90 | 93 |
| 122 | 78 | 87 | 61 | 48 |
| 123 | 95 | 103 | 100 | 105 |
| 124 | 99 | 97 | 104 | 94 |
| 125 | 101 | 102 | 81 | 81 |
| 126 | 102 | 107 | 106 | 106 |
| 127 | 90 | 100 | 94 | 99 |
| 128 | 80 | 75 | 75 | 37 |
| 129 | 95 | 98 | 100 | 99 |
| 130 | 75 | 99 | 84 | 77 |
| 131 | 118 | 142 | 96 | 86 |
| 132 | 94 | 86 | 93 | 62 |

Table 8a: Evaluation of the Questionnaire:

| Lab ID\ Question | Give a reference to your method | Extraction solvent used | Extraction solvent to sample ratio used during extraction (in mL/g)? |
|------------------|--|---|--|
| 101 | SOP | MeOH/H ₂ O=4:1 | 4:1 |
| 102 | Aflaprep | MeOH/H ₂ O=6:4 | 5:1 |
| 103 | EN 14123 | MeOH/H ₂ O=4:1 | 4:1 |
| 104 | EN 14123:2003 | MeOH/H ₂ O=4:1 | 0.167g/mL |
| 105 | J AOAC vol.83/2 N°2, 2000 | MeOH/H ₂ O=4:1 | 4:1 |
| 106 | In house method, LC-MS/MS | AcCN/H ₂ O=4:1 | 4:1 |
| 107 | SR EN 14123/2003 | MeOH/H ₂ O=7:3 | 5:1 |
| 108 | SOP | MeOH | 100mL/28.6g |
| 109 | VICAM Afla Test WB Instruction Manual | MeOH/H ₂ O=7:3 | 5:1 |
| 110 | AOAC 991.31 | MeOH/H ₂ O=7:3 | 5:1 |
| 111 | ISO 16050 | MeOH/H ₂ O=7:3 | 5:1 |
| 112 | Modified EN 14123 | MeOH/H ₂ O=4:1 | 4:1 |
| 113 | PN-EN-14123:2004 | MeOH/H ₂ O=4:1 | 4:1 |
| 114 | in-house developed | AcCN/MeOH=1:1 | 25:3 |
| 115 | EN14123 (2003) | AcCN/MeOH/H ₂ O=1:1:1 | |
| 116 | EN 14123 | 100 ml MeOH/H ₂ O(4/1)+50 ml Hexan | 6:1 |
| 117 | in-house (according EN 14123 | MeOH/H ₂ O=4:1 | 4:1 |
| 118 | EN ISO 17375:2006 | MeOH/H ₂ O=4:1 | 5:1 |
| 119 | Modified AOAC Official Method 991.31 | Acetone/H ₂ O=85/15 | 5:1 |
| 120 | EN 14123 | MeOH/H ₂ O 62.5% | 4:1 |
| 121 | CEN/TC 275 EN 14123:2003 | MeOH/H ₂ O | 4:1 |
| 122 | LVS EN 14123:2003 | MeOH/H ₂ O=4:1 | 4:1 |
| 123 | J.Chromatogr, 1991, 543, 220-225 | MeOH/H ₂ O | 0,25 mL/g |
| 124 | Internal SOP | AcCN/H ₂ O=6:4 | 5:1 |
| 125 | §64-LFGB L 48.00-1 | CHCl ₃ | 5:1 |
| 126 | J. AOAC Int. 1994, 77 (1), 46-53 | Acetone/H ₂ O=85:15 | 6:1 |
| 127 | Instruction manual for the columns MultiSep 226 | AcCN/H ₂ O=4:1 | 5:1 |
| 128 | AflaZon + (Romer lab). | Acetonitrile 84 % in water | 4:1 |
| 129 | CEN standard: prEN 14123 | MeOH/H ₂ O=4:1 | 4:1 |
| 130 | JAOAC Int., 88, 2005, 526 – 535 | MeOH/H ₂ O=6:4 | 5:1 |
| 131 | Method for determination of Aflatoxin B1 and total Aflatoxins in peanut butter, raw peanuts and corn | MeOH/H ₂ O=7:3 | 5:1 |
| 132 | Project SMT-CT96-2045 | MeOH/H ₂ O(4/1)+ Hexan | 6:1 |
| | Senyuva H.Z., Gilbert J., J AOAC Int Vol88, No2, 2005 | H ₂ Oand MeOH | 5:1 |

Table 8b: Evaluation of the Questionnaire:

| Lab ID\ Question | extraction aids added | extraction mode and time | type of clean-up |
|------------------|--|---|---------------------|
| 101 | 0.5g NaCl | shaking 16 hours | IAC |
| 102 | 4g NaCl | Shaking 30 minutes | IAC |
| 103 | 5g NaCl | Blending 3 min | IAC |
| 104 | 2.5g NaCl | Blending 3 min | IAC |
| 105 | 25gNaCl/L | 3min blending + 30min shaking | IAC |
| 106 | No | Shaking for 2 hours | none |
| 107 | 5g NaCl | Blending 3min | IAC |
| 108 | 5g NaCl | Shaking | IAC |
| 109 | 5g NaCl | Blending 5 min | IAC |
| 110 | 5g NaCl to 25g | Blend 2 min | IAC |
| 111 | 5 g of NaCl | Shaking 10 min | IAC |
| 112 | 0,1g NaCl/g | Blending 3 min | IAC |
| 113 | 2,5 g NaCl/25 g | Blending 3 min | IAC |
| 114 | No | Blending:3 min solvent1, 2 min solvent2 | IAC |
| 115 | 2,5 g NaCl to 25 g sample | Shaking 30 min | IAC |
| 116 | 5.0 g NaCl | Blending 2 min | IAC |
| 117 | 1g NaCl/10g sample | blending 3 min | IAC |
| 118 | No | Shaking 1hr | IAC |
| 119 | 10% w/w NaCl to dry sample | Blending 2.5 min | IAC |
| 120 | 0,1g NaCl/g sample | Shaking 30 min | IAC |
| 121 | 5g NaCl and 100 ml n-hexane | Blending 3 min | IAC |
| 122 | No | Shaking 30 min | Romer Labs AflaStar |
| 123 | No | Blending 4 min | IAC |
| 124 | 10g acid washed Celite 545/20g sample/100ml chloroform | Shaking 30 min | IAC |
| 125 | 5g NaCl | blending 3 min | IAC |
| 126 | No | Shaking 1 hr | IAC |
| 127 | No | Shaking 30 min | MultiSep-column |
| 128 | 5g NaCl | Blending 2 min | IAC |
| 129 | 2g NaCl | Blending 1min after addition of H2O and 2min after addition of MeOH | IAC |
| 130 | 5g NaCl | Blending 2 min | IAC |
| 131 | 5g NaCl | Blending 3min | IAC |
| 132 | 2g NaCl | Blending | IAC |

Table 8c: Evaluation of the Questionnaire:

| Lab ID\ Question | Extract evaporated prior injection | derivatisation method applied? If yes, please state the kind of method | any "over-night" stops in the analysis? If yes, please state at what point. |
|------------------|---|--|---|
| 101 | No | PBPB | No |
| 102 | Yes to 1mL | Iodine derivatisation | YES, there were a lot of measurements and the system with the external pump for pumping iodine can not be left alone during night. |
| 103 | Yes | Kobra cell | Yes before injections |
| 104 | No | Kobra cell | Yes, problem with Kobra cell |
| 105 | No | Kobra cell | Yes, After IA-clean-up |
| 106 | No | No | No |
| 107 | No | Kobra cell | No |
| 108 | No | Kobra cell | Yes, after extraction of the samples (28.6g with 100mL Methanol) and filtration of extracts, the extracts were stored in a freezer overnight. |
| 109 | Yes | TFA | No |
| 110 | No | Iodine derivatisation | No |
| 111 | No | Kobra cell | No |
| 112 | No | PBPB | No |
| 113 | No | PBPB | No |
| 114 | Yes | Kobra cell | No |
| 115 | No | PBPB | No |
| 116 | Yes | PBPB | No |
| 117 | No | Kobra cell | Yes, before clean-up |
| 118 | No | Kobra cell | No |
| 119 | No | Iodine derivatisation | No |
| 120 | Yes | TFA | No |
| 121 | No | PBPB | No |
| 122 | Yes | Kobra cell | Yes, Evaporation |
| 123 | No | Kobra cell | No |
| 124 | Yes, under N2 | Kobra cell | No |
| 125 | Yes, under N2 | Kobra cell | No |
| 126 | No | Kobra cell | No |
| 127 | Yes, 2 ml of extract was evaporated, dried and dissolved in 300 µl. | Kobra cell | Yes, after clean-up on MultiSep-column, the 2 ml extract for evaporation was frozen overnight |
| 128 | Yes | PBPB | Stop after elution from IAC |
| 129 | No | PBPB | No |
| 130 | Yes | PBPB | No |
| 131 | No | PBPB | No |
| 132 | No | Kobra cell | No |

Table 8d: Evaluation of the Questionnaire:

| Lab ID\ Question | use acid washed glass ware | protection against daylight | unusual observations |
|------------------|---|------------------------------|--|
| 101 | Yes | Yes | |
| 102 | Yes | Only at measurement by HPLC | No |
| 103 | Yes | Yes | |
| 104 | No | Yes, problem with Kobra cell | No |
| 105 | No | No | Blank material analysed unspiked was shown to contain very low concentrations of Aflatoxin |
| 106 | Yes | Yes | The signal for Aflatoxin B1 of the standard solution (10 times diluted) is out of the calibration range. |
| 107 | Yes | Yes | No |
| 108 | Yes, after extraction of the samples (28.6g with 100mL Methanol) and filtration of extracts, the extracts were stored in a freezer overnight. | Yes | No |
| 109 | Yes | Yes | |
| 110 | No | Yes | No |
| 111 | Yes | Yes | trace of aflatoxin B1 detected in the blank sample Though the requested storage temperature was 4°, the temperature written on containers was -20°, so the samples were stored for some time in the freezer |
| 112 | Yes | Yes | No |
| 113 | Yes | Yes | No |
| 114 | No | Yes | The 'positive' sample was a different colour from the 'blank'. Also it had a different (less 'fresh') odour |
| 115 | Yes | Yes | No |
| 116 | Yes | Yes | No |
| 117 | No | Yes | No |
| 118 | Yes | Yes | Yes, Low recovery for G1 and G2 |
| 119 | Yes | Yes | No |
| 120 | No | No | No |
| 121 | No | Yes-Al foil | No |
| 122 | No | No | problem during extraction, sample material sticks to container |
| 123 | Yes | Yes | AFB1 was outside the calibration range prepared, therefore, further dilutions were made of the standard (x50) and AFB1 concentration was calculated from this |
| 124 | No | Yes | The samples contained a lot of fat, after chloroform extraction it is evaporated and re-dissolved in methanol a striking colour difference between the extract of the blank and the positive sample. |
| 125 | No | Yes | Yes, it was the first time we observed a second fluid layer after extraction and centrifugation. We didn't analysed that solution for aflatoxins! |
| 126 | Yes | Yes | Some sample amount was always left on the bottom of the extraction glassware. |
| 127 | No | Yes | The blank sample for spiking gave peaks at the retention times of toxins. Blank subtraction was made. |
| 128 | No, but 0.1% acetic acid is added to the mobile phase, which is used for solution of standards and samples | Yes | No |
| 129 | No | Yes | No |
| 130 | Yes | Yes | No |
| 131 | Hypochlorite and Acetone/H2O | Yes | The 'positive' sample was a different colour from the 'blank'. |
| 132 | No | Yes | for two blank samples G2 recovery were below 85% |

Table 9a: z-scores for Aflatoxin B1

| LAB ID | Af positive | Test Solution | Af positive TSC | Af positive REC | Af positive CTSREC |
|--------|-------------|---------------|-----------------|-----------------|--------------------|
| 101 | 0.9 | -0.3 | 1.2 | 2.2 | 2.6 |
| 102 | 1.7 | 0.0 | 1.6 | 1.7 | 1.7 |
| 103 | -1.0 | -0.2 | -0.8 | -0.2 | 0.0 |
| 104 | -1.3 | 0.0 | -1.3 | 0.0 | 0.0 |
| 105 | -1.1 | -0.4 | -0.8 | -0.8 | -0.4 |
| 106 | 1.3 | -0.8 | 2.5 | 1.2 | 2.4 |
| 107 | 15.1 | 3.9 | 6.0 | 16.2 | 6.6 |
| 108 | -3.3 | -0.5 | -3.1 | -3.2 | -3.0 |
| 109 | -0.2 | 3.6 | -2.1 | 0.9 | -1.5 |
| 110 | -1.1 | -1.0 | -0.2 | -0.1 | 1.1 |
| 111 | -0.1 | 0.9 | -0.8 | 0.3 | -0.5 |
| 112 | -0.6 | -0.1 | -0.5 | -0.4 | -0.3 |
| 113 | -0.8 | -0.2 | -0.7 | 0.1 | 0.3 |
| 114 | 15.7 | -0.2 | 16.8 | 18.4 | 19.6 |
| 115 | -3.1 | 0.0 | -3.1 | -3.0 | -3.0 |
| 116.1 | -1.1 | -0.3 | -0.9 | -1.0 | -0.8 |
| 116.2 | -1.1 | -0.2 | -1.0 | -1.0 | -0.8 |
| 117 | -2.1 | 0.1 | -2.2 | -1.4 | -1.5 |
| 118 | -0.2 | 0.0 | -0.2 | -0.3 | -0.3 |
| 119 | 0.0 | 0.7 | -0.6 | 0.5 | -0.1 |
| 120 | -3.7 | -0.3 | -3.7 | -3.0 | -3.0 |
| 121 | -0.4 | -0.5 | 0.0 | 0.5 | 1.0 |
| 122 | -0.6 | 0.9 | -1.3 | 0.5 | -0.4 |
| 123 | 0.0 | 0.7 | -0.6 | 0.3 | -0.4 |
| 124 | -0.5 | 0.3 | -0.7 | -0.4 | -0.7 |
| 125 | 1.0 | 1.0 | 0.0 | 1.0 | -0.1 |
| 126 | 0.2 | -0.2 | 0.4 | 0.1 | 0.3 |
| 127 | 0.1 | 0.0 | 0.1 | 0.6 | 0.6 |
| 128 | -1.2 | -0.4 | -0.9 | -0.3 | 0.0 |
| 129 | -0.1 | -0.3 | 0.2 | 0.2 | 0.5 |
| 130 | -1.1 | 3.7 | -2.6 | 0.1 | -2.0 |
| 131 | 13.1 | 6.2 | 2.9 | 10.4 | 1.8 |
| 132 | -0.8 | -0.7 | -0.2 | -0.6 | 0.1 |

Af positive = Score rating on the neat result for the Af positive sample, prior any correction; Test Solution = Score rating for the result for the test solution; Af positive TSC = Score rating on the result of the Af positive sample, after correction with the test solution result; Af positive REC = Score rating on the result of the Af positive sample, after correction for recovery; Af positive CRC = Score rating on the result of the Af positive sample, after correction with the test solution result and recovery.

Table 9b: z-scores for Aflatoxin B2

| LAB ID | Af positive | Test Solution | Af positive TSC | Af positive REC | Af positive CTSREC |
|--------|-------------|---------------|-----------------|-----------------|--------------------|
| 101 | 0.1 | 0.0 | 0.2 | 1.1 | 1.2 |
| 102 | 3.4 | 1.9 | 1.0 | 3.4 | 1.0 |
| 103 | -0.7 | -0.1 | -0.6 | 0.2 | 0.3 |
| 104 | -0.6 | 0.0 | -0.6 | 0.6 | 0.6 |
| 105 | -0.7 | 0.2 | -0.9 | -0.1 | -0.3 |
| 106 | -4.5 | -0.4 | -4.5 | -4.5 | -4.5 |
| 107 | 20.2 | 5.2 | 6.9 | 23.7 | 8.6 |
| 108 | -3.2 | 0.7 | -3.4 | -3.1 | -3.3 |
| 109 | -1.2 | 27.3 | -4.1 | 0.4 | -3.8 |
| 110 | -1.0 | -0.7 | -0.4 | 0.0 | 0.8 |
| 111 | -1.9 | -0.4 | -1.6 | 0.2 | 0.7 |
| 112 | -0.3 | 0.2 | -0.5 | 0.2 | -0.1 |
| 113 | -0.3 | -0.1 | -0.2 | 0.6 | 0.7 |
| 114 | 17.8 | 0.1 | 17.4 | 20.6 | 20.1 |
| 115 | -3.4 | 0.5 | -3.6 | -3.3 | -3.5 |
| 116.1 | -0.8 | 0.0 | -0.8 | -0.7 | -0.8 |
| 116.2 | -1.1 | 0.0 | -1.1 | -0.9 | -1.0 |
| 117 | -2.4 | 0.1 | -2.5 | -1.6 | -1.6 |
| 118 | 0.3 | 0.1 | 0.2 | -0.1 | -0.2 |
| 119 | 0.3 | 0.8 | -0.5 | 0.6 | -0.2 |
| 120 | -3.6 | -0.3 | -3.5 | -2.9 | -2.8 |
| 121 | 0.5 | -0.2 | 0.7 | 0.9 | 1.1 |
| 122 | -0.3 | 0.5 | -0.7 | 0.4 | -0.1 |
| 123 | 0.9 | 1.1 | -0.2 | 0.8 | -0.3 |
| 124 | -0.3 | 0.3 | -0.5 | -0.1 | -0.4 |
| 125 | 1.5 | 1.0 | 0.4 | 1.4 | 0.3 |
| 126 | 0.3 | -0.1 | 0.4 | 0.0 | 0.1 |
| 127 | 1.0 | 0.3 | 0.6 | 1.0 | 0.6 |
| 128 | -1.3 | -0.7 | -0.8 | -0.3 | 0.4 |
| 129 | -0.3 | -0.7 | 0.5 | -0.2 | 0.6 |
| 130 | 0.3 | 7.3 | -2.7 | 0.3 | -2.7 |
| 131 | 8.0 | 3.6 | 2.5 | 4.3 | 0.4 |
| 132 | -0.2 | -0.3 | 0.1 | 0.5 | 0.9 |

Af positive = Score rating on the neat result for the Af positive sample, prior any correction; Test Solution = Score rating for the result for the test solution; Af positive TSC = Score rating on the result of the Af positive sample, after correction with the test solution result; Af positive REC = Score rating on the result of the Af positive sample, after correction for recovery; Af positive CRC = Score rating on the result of the Af positive sample, after correction with the test solution result and recovery.

Table 9c: z-scores for Aflatoxin G1

| LAB ID | Af positive | Test Solution | Af positive TSC | Af positive REC | Af positive CTSREC |
|--------|-------------|---------------|-----------------|-----------------|--------------------|
| 101 | -1.5 | 0.3 | -1.7 | -0.6 | -0.9 |
| 102 | -0.4 | 2.0 | -1.7 | 1.8 | -0.1 |
| 103 | -2.3 | -0.3 | -2.1 | -0.2 | 0.0 |
| 104 | -1.3 | 0.1 | -1.4 | -0.4 | -0.5 |
| 105 | -1.2 | 0.2 | -1.3 | -0.7 | -0.9 |
| 106 | 1.1 | -1.9 | 5.2 | 0.7 | 4.6 |
| 107 | 15.1 | 5.2 | 4.6 | 15.9 | 5.0 |
| 108 | -2.8 | 0.2 | -2.9 | -2.6 | -2.7 |
| 109 | 0.5 | 14.0 | -3.3 | 1.6 | -3.0 |
| 110 | -0.9 | 0.0 | -0.9 | 0.3 | 0.3 |
| 111 | -0.3 | 1.0 | -1.1 | 3.1 | 1.7 |
| 112 | -1.2 | 0.2 | -1.4 | -0.6 | -0.8 |
| 113 | -1.8 | -1.3 | -0.8 | -1.1 | 0.1 |
| 114 | 12.7 | 0.1 | 12.4 | 15.7 | 15.3 |
| 115 | -2.8 | 0.3 | -2.9 | -2.8 | -2.9 |
| 116.1 | -1.9 | -0.4 | -1.7 | -1.5 | -1.2 |
| 116.2 | -2.1 | -0.7 | -1.7 | -1.7 | -1.2 |
| 117 | -2.4 | -0.2 | -2.3 | -1.6 | -1.5 |
| 118 | -2.8 | 0.2 | -2.9 | -0.8 | -1.0 |
| 119 | 0.0 | 1.0 | -0.8 | 0.5 | -0.4 |
| 120 | -3.7 | 0.1 | -3.7 | -2.8 | -2.8 |
| 121 | -0.2 | 0.2 | -0.3 | 0.3 | 0.1 |
| 122 | -1.6 | 1.0 | -2.1 | 0.3 | -0.6 |
| 123 | 0.1 | 0.5 | -0.4 | 0.1 | -0.4 |
| 124 | -0.7 | 0.2 | -0.9 | -0.8 | -1.0 |
| 125 | 0.5 | 1.6 | -0.9 | 1.6 | 0.0 |
| 126 | 0.2 | -0.2 | 0.3 | -0.1 | 0.1 |
| 127 | 0.1 | 0.4 | -0.2 | 0.5 | 0.1 |
| 128 | -1.7 | -0.5 | -1.4 | -0.7 | -0.3 |
| 129 | 1.4 | 1.0 | 0.2 | 1.3 | 0.2 |
| 130 | -0.9 | 6.4 | -3.0 | -0.2 | -2.7 |
| 131 | 11.7 | 7.8 | 1.4 | 12.4 | 1.7 |
| 132 | 0.5 | -0.4 | 0.9 | 0.9 | 1.3 |

Af positive = Score rating on the neat result for the Af positive sample, prior any correction; Test Solution = Score rating for the result for the test solution; Af positive TSC = Score rating on the result of the Af positive sample, after correction with the test solution result; Af positive REC = Score rating on the result of the Af positive sample, after correction for recovery; Af positive CRC = Score rating on the result of the Af positive sample, after correction with the test solution result and recovery.

Table 9d: z-scores for Aflatoxin G2

| LAB ID | Af positive | Test Solution | Af positive TSC | Af positive REC | Af positive CTSREC |
|--------|-------------|---------------|-----------------|-----------------|--------------------|
| 101 | -1.7 | 0.1 | -1.7 | -0.5 | -0.5 |
| 102 | 1.2 | 0.7 | 0.4 | 6.3 | 4.8 |
| 103 | -2.9 | 0.1 | -3.0 | 2.4 | 2.2 |
| 104 | -1.9 | 0.0 | -1.9 | -1.0 | -1.1 |
| 105 | -2.3 | 0.1 | -2.4 | 1.5 | 1.4 |
| 106 | -4.5 | -1.9 | -4.5 | -4.5 | -4.5 |
| 107 | 10.4 | 5.1 | 2.5 | 35.1 | 14.2 |
| 108 | -2.9 | 0.4 | -3.1 | -2.7 | -2.9 |
| 109 | -2.0 | 34.3 | -4.2 | -0.9 | -4.1 |
| 110 | -4.5 | -0.6 | -4.5 | -4.5 | -4.5 |
| 111 | -2.6 | 0.1 | -2.7 | -1.3 | -1.4 |
| 112 | -1.5 | 0.4 | -1.7 | -0.9 | -1.2 |
| 113 | -0.3 | 0.6 | -0.9 | 1.1 | 0.4 |
| 114 | 13.3 | 0.2 | 12.6 | 15.9 | 15.0 |
| 115 | -3.1 | 0.6 | -3.2 | -2.3 | -2.5 |
| 116.1 | -1.8 | 0.2 | -1.9 | -0.6 | -0.7 |
| 116.2 | -1.6 | 0.4 | -1.8 | -0.6 | -0.9 |
| 117 | -2.0 | 0.2 | -2.1 | -1.0 | -1.1 |
| 118 | -3.1 | 0.9 | -3.4 | 1.3 | 0.3 |
| 119 | -1.6 | 0.1 | -1.7 | -1.2 | -1.2 |
| 120 | -3.6 | 0.6 | -3.7 | -2.7 | -2.9 |
| 121 | -0.6 | 1.1 | -1.4 | -0.3 | -1.1 |
| 122 | -2.5 | 0.8 | -2.8 | -0.3 | -0.9 |
| 123 | -0.2 | 0.7 | -0.8 | -0.4 | -1.0 |
| 124 | -0.9 | 1.3 | -1.7 | -0.7 | -1.6 |
| 125 | 1.1 | 2.8 | -1.1 | 2.5 | -0.2 |
| 126 | -0.5 | 0.1 | -0.6 | -0.7 | -0.8 |
| 127 | 0.7 | 0.5 | 0.1 | 0.7 | 0.2 |
| 128 | -2.2 | 0.1 | -2.2 | 1.8 | 1.7 |
| 129 | -0.8 | -0.1 | -0.7 | -0.8 | -0.7 |
| 130 | -3.3 | 8.7 | -4.1 | -3.0 | -4.0 |
| 131 | 4.1 | 7.4 | -1.3 | 5.5 | -0.7 |
| 132 | 0.1 | 0.8 | -0.6 | 3.0 | 1.8 |

Af positive = Score rating on the neat result for the Af positive sample, prior any correction; Test Solution = Score rating for the result for the test solution; Af positive TSC = Score rating on the result of the Af positive sample, after correction with the test solution result; Af positive REC = Score rating on the result of the Af positive sample, after correction for recovery; Af positive CRC = Score rating on the result of the Af positive sample, after correction with the test solution result and recovery.

Table 9e: z-scores for Total Aflatoxin

| LAB ID | Af positive | Test Solution | Af positive TSC | Af positive REC | Af positive CTSREC |
|--------|-------------|---------------|-----------------|-----------------|--------------------|
| 101 | -0.1 | - | 0.0 | 1.1 | 1.2 |
| 102 | 1.3 | - | 0.6 | 2.4 | 1.4 |
| 103 | -1.5 | - | -1.3 | 0.1 | 0.3 |
| 104 | -1.3 | - | -1.3 | -0.1 | -0.2 |
| 105 | -1.2 | - | -1.1 | -0.5 | -0.4 |
| 106 | -0.1 | - | 1.6 | -0.2 | 1.4 |
| 107 | 15.4 | - | 5.4 | 18.9 | 7.1 |
| 108 | -3.1 | - | -3.1 | -3.0 | -3.0 |
| 109 | -0.3 | - | -2.9 | 0.8 | -2.5 |
| 110 | -1.4 | - | -0.8 | -0.4 | 0.3 |
| 111 | -0.6 | - | -1.2 | 0.9 | 0.2 |
| 112 | -0.8 | - | -0.8 | -0.4 | -0.5 |
| 113 | -1.0 | - | -0.7 | -0.1 | 0.3 |
| 114 | 15.0 | - | 15.4 | 17.7 | 18.2 |
| 115 | -3.1 | - | -3.1 | -2.9 | -3.0 |
| 116.1 | -1.4 | - | -1.2 | -1.0 | -0.9 |
| 116.2 | -1.3 | - | -1.3 | -1.1 | -0.9 |
| 117 | -2.2 | - | -2.2 | -1.5 | -1.5 |
| 118 | -1.1 | - | -1.1 | -0.2 | -0.4 |
| 119 | -0.1 | - | -0.7 | 0.4 | -0.3 |
| 120 | -3.7 | - | -3.7 | -2.9 | -2.9 |
| 121 | -0.3 | - | -0.1 | 0.4 | 0.6 |
| 122 | -1.0 | - | -1.6 | 0.4 | -0.4 |
| 123 | 0.1 | - | -0.5 | 0.2 | -0.4 |
| 124 | -0.5 | - | -0.8 | -0.5 | -0.8 |
| 125 | 1.0 | - | -0.3 | 1.3 | 0.0 |
| 126 | 0.2 | - | 0.3 | 0.0 | 0.1 |
| 127 | 0.3 | - | 0.1 | 0.6 | 0.4 |
| 128 | -1.4 | - | -1.1 | -0.2 | 0.2 |
| 129 | 0.2 | - | 0.2 | 0.3 | 0.3 |
| 130 | -1.0 | - | -2.9 | -0.2 | -2.5 |
| 131 | 11.2 | - | 2.1 | 9.6 | 1.3 |
| 132 | -0.3 | - | 0.1 | 0.3 | 0.7 |

Af positive = Score rating on the neat result for the Af positive sample, prior any correction; Test Solution = Score rating for the result for the test solution; Af positive TSC = Score rating on the result of the Af positive sample, after correction with the test solution result; Af positive REC = Score rating on the result of the Af positive sample, after correction for recovery; Af positive CRC = Score rating on the result of the Af positive sample, after correction with the test solution result and recovery.

EUR 23261 EN– Joint Research Centre – Institute for Reference Materials and Measurements

Title: Report on the 2007 Proficiency Test of the Community Reference Laboratory Network – Determination of Aflatoxins in a Peanut Product and a Test Solution

Author(s): STROKA, Joerg, BREIDBACH, Andreas, DONCHEVA TSANEVA, Ivanka, AMBROSIO, Massimo
Luxembourg: Office for Official Publications of the European Communities

2008– 31 pp. – 21 x 29.7 cm

EUR – Scientific and Technical Research series – ISSN 1018-5593

ISBN 978-91-79-08394-5

DOI 10.2787/28654

Abstract

A proficiency test was conducted with 31 European National Reference Laboratories (NRLs) for mycotoxins and one Laboratory from Turkey. Test materials were a mixed aflatoxin (Af) solution in acetonitrile and two candidate Certified Reference Materials ("aflatoxin positive" and "blank") that have not yet been released. Laboratories determined the aflatoxin content by reverse-phase high-performance liquid-chromatography (RP-HPLC) with either fluorescence or mass-selective detection against their own standard solutions as reference.

Applying the modified Horwitz equation according to Thompson as a basis for the target standard deviation (22% in the case of this proficiency test), 26 out of 32 NRLs achieved z-scores of less than 2 and 17 NRLs reported values within the uncertainty range for both aflatoxin B1 and total aflatoxins in the candidate CRM after correction for recovery in both cases.

How to obtain EU publications

Our priced publications are available from EU Bookshop (<http://bookshop.europa.eu>), where you can place an order with the sales agent of your choice.

The Publications Office has a worldwide network of sales agents. You can obtain their contact details by sending a fax to (352) 29 29-42758.

The mission of the JRC is to provide customer-driven scientific and technical support for the conception, development, implementation and monitoring of EU policies. As a service of the European Commission, the JRC functions as a reference centre of science and technology for the Union. Close to the policy-making process, it serves the common interest of the Member States, while being independent of special interests, whether private or national.



ISBN 978-92-79-08394-5

